



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

Handwritten notes: TSG, 4.11.2012

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/29/12 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 3/1/12 and 3/29/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Handwritten signature: Sarah Cione

Sarah L. Cione


Handwritten signature: 4/9/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th day of April, 2012.

Handwritten signature: Adina Escarcega

Notary Public

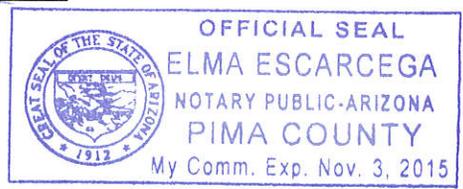

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

Bg
4.11.2012

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/29/12 TIME 1439

INTOXILYZER SERIAL # 80-002262 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
03/29/2012

14:41:16

2 of 7_{sc}

B4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:42:00
Diagnostic Test	Pass	14:42:27
Air Blank	0.000	14:42:55
0.100 Cal Check	0.100	14:43:13
Air Blank	0.000	14:43:43
Subject Test	0.000	14:44:10
Air Blank	0.000	14:44:38
Five Minute Wait	ABT*	
Air Blank	0.000	14:45:17

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
03/29/2012

14:46:23

PG

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:46:56
Air Blank	IPS**	14:47:04

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002262
03/29/2012 14:54:08

TS4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:54:50
Diagnostic Test	Pass	14:55:17
Air Blank	0.000	14:55:45
0.100 Cal Check	0.100	14:56:03
Air Blank	0.000	14:56:33
Subject Test	0.139*	14:59:41
Air Blank	0.000	15:00:13
Air Blank	0.000	15:00:40
Subject Test	XXX**	15:01:08
Air Blank	0.000	15:01:42

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7_{sc}


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
03/29/2012 15:01:59

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:02:26
Cal Check	0.099	15:02:45
Air Blank	0.000	15:03:14

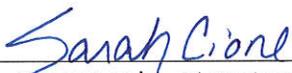

Operator's Signature

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
4.11.2012

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/29/12 TIME 1439
Intoxilyzer Serial # 80-002262 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002262
03/29/2012 14:39:12

7 of 7_{sc}

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:39:39
Cal Check	0.099	14:39:58
Air Blank	0.000	14:40:27

Sarah Cione

Operator's Signature