



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
8-2-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/31/12 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of - and - the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on 7/31/12.

Sarah Cione

Sarah L. Cione


8/2/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 2nd, day of August, 2012.

Elma Escarcega

Notary Public

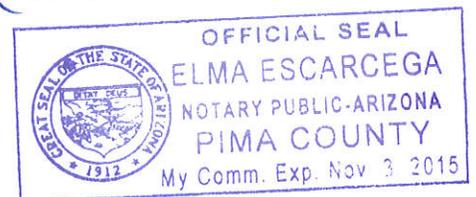

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 12 sc

INTOXILYZER SERVICE RECORD

TS4  
8-2-2012

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002262

Date Returned to Service: 7/31/12

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
8-2-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/31/12 TIME 1431

INTOXILYZER SERIAL # 80-002262 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
07/31/2012

14:33:23

3 of 12 sc  
Bh

---

Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:34:07
Diagnostic Test	Pass	14:34:42
Air Blank	0.000	14:35:09
0.100 Cal Check	0.100	14:35:30
Air Blank	0.000	14:35:59
Subject Test	0.000	14:36:13
Air Blank	0.000	14:36:40
Wait	ABT*	
Air Blank	0.000	14:37:17

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

4 of 12 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
07/31/2012 14:37:42

---

Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:38:26
Diagnostic Test	Pass	14:39:01
Air Blank	0.000	14:39:28
0.100 Cal Check	0.100	14:39:49
Air Blank	0.000	14:40:18
Subject Test	RFI*	14:40:34
Air Blank	IPS**	14:40:47

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

5 of 12 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
07/31/2012 14:40:51

---

Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:41:37
Diagnostic Test	Pass	14:42:12
Air Blank	0.000	14:42:39
0.100 Cal Check	0.101	14:43:00
Air Blank	0.000	14:43:29
Subject Test	0.077*	14:46:37
Air Blank	0.000	14:47:09
Air Blank	0.000	14:47:36
Subject Test	XXX**	14:48:06
Air Blank	0.000	14:48:39

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

6 of 12 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
07/31/2012 14:49:00

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:49:28
0.100 Cal Check	0.100	14:49:49
Air Blank	0.000	14:50:19

Preventative Maintenance NOT  
Performed -- Timer NOT Reset

EXHIBIT G-2

7 of 12 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

TS4  
8-2-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 7/31/12 TIME 1431  
Intoxilyzer Serial # 80-002262 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

8 of 12 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
07/31/2012 14:31:43

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:32:14
0.100 Cal Check	0.100	14:32:35
Air Blank	0.000	14:33:05

Preventative Maintenance  
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 12 sc  
TS4  
8-2-2012

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 7/31/12  
 QAS: Sarah Ciont  
 Intoxilyzer Serial Number: 80-002262  
 New Lot Number: 09312100A1  
 New Expiration Date: 6/1/2014

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
 Select "D" for dry gas and enter the following information:  
     Enter Target Value: 0.100  
     Enter Control Lot #: From the tank label  
     Enter Exp Date: From the tank label (mm/dd/yyyy)  
     Edit info Y/N? Y or N  
     Update Std Info? Y ENTER  
     "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove ~~and sign~~ printed record. Attach the record to the completed checklist.  
     sc

SIGNATURE Sarah Ciont

10 of 12 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
07/31/2012 14:17:31

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:18:01
0.100 Cal Check	0.100	14:18:22
Air Blank	0.000	14:18:52

Preventative Maintenance  
Performed -- Timer Reset



12 of 12 sc  
134

## Certificate of Calibration

This is to certify the calibration of **Intoxilyzer**® serial number 80-002262, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to conform to the National Highway Traffic Safety Administration (NHTSA) Standard for Devices to Measure Breath Alcohol (Federal Register, Vol.58, No.179, pp 48705-48710, Sept. 17, 1993) for accuracy and precision. Reference materials are traceable through the National Institute of Standards and Technology (NIST) to the International System of Units (SI).

Date July 23, 2012 Signed D. R. R. R.  
Technician

**EMI** INC.

316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650517 Rev.A

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.