



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 4.17.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/28/13 Cal 1 0.103 AC Cal 2 0.102 AC

Between the dates of 2/28/13 and 3/28/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.


Sarah Cione
 Sarah L. Cione

4/8/13
 ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 8th, day of April, 2013.

Elma Escarcega
 Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

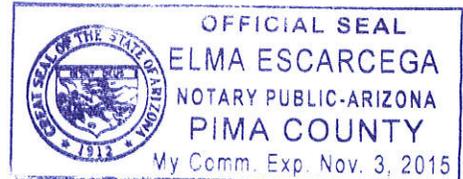

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

TSK
4.17.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/28/13 TIME 759

INTOXILYZER SERIAL # 80-002262 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7sc
TBH

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/28/2013 08:01:35

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:02:21
Diagnostic Test	Pass	08:02:56
Air Blank	0.000	08:03:23
0.100 Cal Check	0.103	08:03:44
Air Blank	0.000	08:04:14
Subject Test	0.000	08:05:07
Air Blank	0.000	08:05:34
Wait	ABT*	
Air Blank	0.000	08:06:14

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_x

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/28/2013

08:12:37

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:13:40
Diagnostic Test	Pass	08:14:15
Air Blank	0.000	08:14:43
0.100 Cal Check	0.102	08:15:03
Air Blank	0.000	08:15:33
Subject Test	RFI*	08:16:07
Air Blank	IPS**	08:16:15

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 SC
D4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/28/2013 08:17:20

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:18:15
Diagnostic Test	Pass	08:18:50
Air Blank	0.000	08:19:17
0.100 Cal Check	0.103	08:19:38
Air Blank	0.000	08:20:07
Subject Test	0.185*	08:23:15
Air Blank	0.000	08:23:47
Air Blank	0.000	08:24:14
Subject Test	XXX**	08:24:39
Air Blank	0.000	08:25:12

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7_{sc}
B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/28/2013 08:26:57

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:27:26
0.100 Cal Check	0.102	08:27:46
Air Blank	0.000	08:28:16

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7
734
4.17.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/28/13 TIME 759
Intoxilyzer Serial # 80-002262 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/28/2013

07:59:16

7 of 7 SC
T84

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:59:44
0.100 Cal Check	0.103	08:00:05
Air Blank	0.000	08:00:35

Preventative Maintenance
Performed -- Timer Reset