



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS9  
2-11-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/30/14 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 1/3/14 and 1/30/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/11/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11<sup>th</sup>, day of February, 2014.

Elma Escarcega

Notary Public

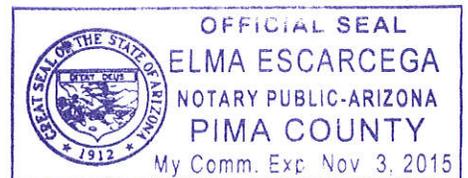

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 SC  
TS4  
2-11-2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/30/14 TIME 932

INTOXILYZER SERIAL # 80-002262 LOCATION DUI HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
01/30/2014

09:34:24

2 of 7  
sc  
T84

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:35:23
Diagnostic Test	Pass	09:35:58
Air Blank	0.000	09:36:26
0.100 Cal Check	0.099	09:36:46
Air Blank	0.000	09:37:16
Subject Test	0.000	09:37:41
Air Blank	0.000	09:38:08
Wait	ABT*	
Air Blank	0.000	09:39:02

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
01/30/2014

09:41:05

3 of 7<sub>sc</sub>  
TD4

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:42:38
Diagnostic Test	Pass	09:43:13
Air Blank	0.000	09:43:40
0.100 Cal Check	0.099	09:44:01
Air Blank	0.000	09:44:30
Subject Test	RFI*	09:44:50
Air Blank	IPS**	09:44:57

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
01/30/2014

09:45:12

4 of 7  
sc  
TS4

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:46:23
Diagnostic Test	Pass	09:46:59
Air Blank	0.000	09:47:26
0.100 Cal Check	0.099	09:47:47
Air Blank	0.000	09:48:16
Subject Test	0.169*	09:51:24
Air Blank	0.000	09:51:56
Air Blank	0.000	09:52:23
Subject Test	XXX**	09:52:49
Air Blank	0.000	09:53:22

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
01/30/2014 09:53:44

5 of 7  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:54:12
0.100 Cal Check	0.098	09:54:33
Air Blank	0.000	09:55:02

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
sc  
TSh  
2.11.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 1/30/14 TIME 932  
Intoxilyzer Serial # 80-002262 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002262

Core Version: 8105.48

01/30/2014

09:32:14

*7 of 7 sc*  
*TS4*

-----  
Standard Lot#: 30412100A1

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Operator: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:32:42
0.100 Cal Check	0.099	09:33:03
Air Blank	0.000	09:33:33

Preventative Maintenance  
Performed -- Timer Reset