



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

SC  
 12/10/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11-6-2014 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 10-9-2014 and 11-6-2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos

12-4-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4<sup>th</sup>, day of December, 2014.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

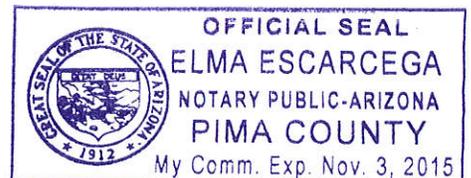


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
TS4  
SC  
12/10/14

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 11-6-2014 TIME 1411

INTOXILYZER SERIAL # 80-002262 LOCATION ODW-DUI-HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check. *off by 1 hr for daylight savings time. Reset clock.*  
(✓) 2. Date check. *TS4*

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
11/06/2014

14:14:56

2 of 7  
TS4

SC

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:15:48
Diagnostic Test	Pass	14:16:23
Air Blank	0.000	14:16:50
0.100 Cal Check	0.100	14:17:11
Air Blank	0.000	14:17:41
Subject Test	0.000	14:18:06
Air Blank	0.000	14:18:33
Wait	ABT*	
Air Blank	0.000	14:19:07

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
11/06/2014

14:26:16

3 of 7  
TS4

50

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:27:03
Diagnostic Test	Pass	14:27:38
Air Blank	0.000	14:28:05
0.100 Cal Check	0.100	14:28:26
Air Blank	0.000	14:28:55
Subject Test	0.077*	14:32:03
Air Blank	0.000	14:32:35
Air Blank	0.000	14:33:02
Subject Test	XXX**	14:33:28
Air Blank	0.000	14:34:02

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
11/06/2014

14:22:23

4 of 7  
TS4

SC

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:23:47
Diagnostic Test	Pass	14:24:22
Air Blank	0.000	14:24:49
0.100 Cal Check	RFI*	14:25:03
Air Blank	IPS**	14:25:15

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
11/06/2014 14:35:34

Sof7  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:36:03
0.100 Cal Check	0.099	14:36:24
Air Blank	0.000	14:36:54

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6087  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
12/10/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 11-6-2014 TIME 1411  
Intoxilyzer Serial # 80-00 2262 Location ODW-DWI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
11/06/2014 14:11:37

7of7  
T54  
SC

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:12:07
0.100 Cal Check	0.100	14:12:28
Air Blank	0.000	14:12:57

Preventative Maintenance  
Performed -- Timer Reset