



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
12/19/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-4-2014 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 11-6-2014 and 12-4-2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Terry S. Gallegos

Terry S. Gallegos


12-17-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 17th, day of December, 2014.

Elma Escarcega

Notary Public

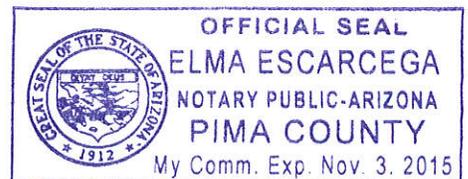

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10 of 9
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
12/19/14

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 12-4-2014 TIME 1457

INTOXILYZER SERIAL # 80-002262 LOCATION ODW-DUI-HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
12/04/2014

15:03:37

2 of 9
TS4

3C

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:04:28
Diagnostic Test	Pass	15:05:03
Air Blank	0.000	15:05:30
0.100 Cal Check	0.101	15:05:51
Air Blank	0.000	15:06:21
Subject Test	0.000	15:06:46
Air Blank	0.000	15:07:13
Wait	ABT*	
Air Blank	0.000	15:07:49

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
12/04/2014

15:08:54

30f9
D4

5C

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:09:47
Diagnostic Test	Pass	15:10:22
Air Blank	0.000	15:10:50
0.100 Cal Check	RFI*	15:11:06
Air Blank	IPS**	15:11:18

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
12/04/2014

15:14:59

4089
754

5C

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:16:21
Diagnostic Test	Pass	15:16:56
Air Blank	0.000	15:17:24
0.100 Cal Check	0.101	15:17:45
Air Blank	0.000	15:18:14
Subject Test	0.047*	15:21:22
Air Blank	0.000	15:21:53
Air Blank	0.000	15:22:20
Subject Test	XXX**	15:23:02
Air Blank	0.000	15:23:37

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
12/04/2014 15:24:05

Sof 9
TS4
SC

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:24:34
0.100 Cal Check	0.101	15:24:54
Air Blank	0.000	15:25:24

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6069
B4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
12/19/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12-4-2014 TIME 1457
Intoxilyzer Serial # 80-002262 Location ODW-DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
12/04/2014 14:57:43

70f9
Bh

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

5c

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:58:12
0.100 Cal Check	0.101	14:58:33
Air Blank	0.000	14:59:02

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

80f9
TS4

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

sc
12/19/14

Date: 12.4.2014
 QAS: Terry Gallegos 28847
 Intoxilyzer Serial Number: 80-002262
 New Lot Number: 11514100A2
 New Expiration Date: 5.1.2016

- (✓) 1. Intoxilyzer 8000 display reads "~~PUSH BUTTON TO START~~" "Contact QAS - Std Change Required" TS4 (expired)
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
12/04/2014 14:49:09

9 of 9
Bh

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

5C

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:49:45
0.100 Cal Check	0.101	14:50:05
Air Blank	0.000	14:50:35

Preventative Maintenance
Performed -- Timer Reset