



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

Alb
 10/23/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10.8.2015 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 9.10.2015 and 10.8.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


10.16.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 16th, day of October, 2015.

Elma Escarcega

Notary Public

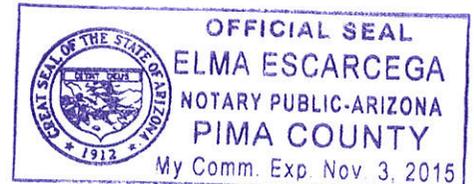

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
B4
Mo
10/23/15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallego 28847 AGENCY Tucson Police Department

DATE 10-8-2015 TIME 1109

INTOXILYZER SERIAL # 80-002262 LOCATION DUI-HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallego

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
10/08/2015 11:23:59

2 of 7
TS4
JL

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:24:52
Diagnostic Test	Pass	11:25:27
Air Blank	0.000	11:25:54
0.100 Cal Check	0.099	11:26:15
Air Blank	0.000	11:26:45
Subject Test	0.000	11:27:06
Air Blank	0.000	11:27:34
Wait	ABT*	
Air Blank	0.000	11:29:07

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
10/08/2015

11:29:13

3 of 7
TS4
MS

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:30:17
Diagnostic Test	Pass	11:30:52
Air Blank	0.000	11:31:19
0.100 Cal Check	RFI*	11:31:33
Air Blank	IPS**	11:31:50

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
10/08/2015

11:32:51

6 of 7
TSY
ADD

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:33:43
Diagnostic Test	Pass	11:34:18
Air Blank	0.000	11:34:45
0.100 Cal Check	0.099	11:35:06
Air Blank	0.000	11:35:35
Subject Test	0.041*	11:38:43
Air Blank	0.000	11:39:14
Air Blank	0.000	11:39:42
Subject Test	XXX**	11:39:58
Air Blank	0.000	11:40:32

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
10/08/2015 11:41:10

50f7
TS4
JLD

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:41:39
0.100 Cal Check	0.099	11:42:00
Air Blank	0.000	11:42:29

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TS4
All
10/23/15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 10.8.2015 TIME 1109
Intoxilyzer Serial # 80-002262 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7077
TS4
dd

INTOXILYZER 8000
Instrument Initialization
10:27:14 10/08/2015

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
10/08/2015 11:09:52

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:10:22
0.100 Cal Check	0.099	11:10:42
Air Blank	0.000	11:11:12

Preventative Maintenance
Performed -- Timer Reset