



WLL  
2/11/16

# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-00 2262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-28-2016 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 12-30-2015 and 1-28-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.



Terry S. Gallegos

2-9-2016


ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 9<sup>th</sup>, day of February, 2016.



Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

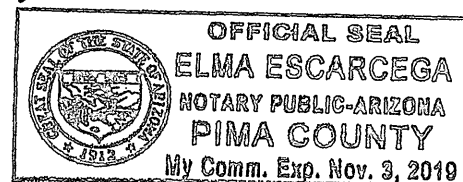

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

10 of 7  
TSG  
WLL  
2/11/16

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 1-28-2016 TIME 1641  
INTOXILYZER SERIAL # 80-002262 LOCATION DUI-HQ

( ☒ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ☒ ) 1. Clock time check.  
( ☒ ) 2. Date check.

OPERATIONAL TESTS

- ( ☒ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ☒ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
( ☒ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ☒ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002262

Core Version: 8105.48

01/28/2016

16:25:14

205.7

B4

WLV

Standard Lot#: 11514100A2

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:26:02
Diagnostic Test	Pass	16:26:37
Air Blank	0.000	16:27:04
0.100 Cal Check	0.101	16:27:25
Air Blank	0.000	16:27:55
Subject Test	0.000	16:28:41
Air Blank	0.000	16:29:08
Wait	ABT*	
Air Blank	0.000	16:29:42

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
01/28/2016

16:30:08

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:31:01
Diagnostic Test	Pass	16:31:36
Air Blank	0.000	16:32:03
0.100 Cal Check	RFI*	16:32:19
Air Blank	IPS**	16:32:32

\*RFI Detect

\*\*Improper Sample

3087  
B4  
WVV

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002262

Core Version: 8105.48

01/28/2016

16:33:03

4 of 7  
B4  
WLL

Standard Lot#: 11514100A2

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:33:51
Diagnostic Test	Pass	16:34:26
Air Blank	0.000	16:34:53
0.100 Cal Check	0.101	16:35:14
Air Blank	0.000	16:35:43
Subject Test	0.111*	16:38:51
Air Blank	0.000	16:39:23
Air Blank	0.000	16:39:51
Subject Test	XXX**	16:40:06
Air Blank	0.000	16:40:39

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
01/28/2016 16:23:25

50f 57  
B4 B4  
WLL

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:23:54
0.100 Cal Check	0.101	16:24:15
Air Blank	0.000	16:24:44

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28842 AGENCY Tucson Police Department  
DATE 1.28.2016 TIME 1641  
Intoxilyzer Serial # 80-002262 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

60f 7  
TS4

WLL  
2/11/16

70f7  
D4  
WLL

INTOXILYZER 8100  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
01/28/2015 16:41:05

Standard Lot#: 1151410042  
Last Changed By: T. GALLEGOS #28847

OPS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Test	g/21.0L	Time
Air Blank	0.000	16:41:33
0.100 Cal Check	0.101	16:41:54
Air Blank	0.000	16:42:24

Preventative Maintenance  
Performed -- Timer Reset