



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
2.11.2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/4/2016 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 1/28/2016 and 2/4/2016  
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis

2/10/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 10<sup>th</sup>, day of February, 2016.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

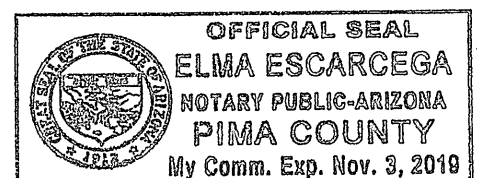


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
Bk 2.11.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 2/4/2016 TIME 09:54  
INTOXILYZER SERIAL # 80-002262 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

William Loomis

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002262

Core Version: 8105.48

02/04/2016

09:56:29

Standard Lot#: 11514100A2

Last Changed By: T. GALLEGOS

#28847

QAS: W. LOOMIS

#53545

TUCSON PD

Operator: W. LOOMIS

#53545

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M

Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:57:20
Diagnostic Test	Pass	09:57:55
Air Blank	0.000	09:58:22
0.100 Cal Check	0.102	09:58:43
Air Blank	0.000	09:59:13
Subject Test	0.000	09:59:30
Air Blank	0.000	09:59:57
Wait	ABT*	
Air Blank	0.000	10:00:32

\*Sequence Aborted

2 of 7  
WR  
134

## Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002262

Core Version: 8105.48

02/04/2016

10:00:35

Standard Lot#: 11514100A2

Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545

TUCSON PD

Operator: W. LOOMIS #53545

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M

Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:01:31
Diagnostic Test	Pass	10:02:06
Air Blank	0.000	10:02:34
0.100 Cal Check	RFI*	10:02:48
Air Blank	IPS**	10:02:57

\*RFI Detect

\*\*Improper Sample

3 of 7  
WLL  
Bh

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-002262

Core Version: 8105.48

02/04/2016

10:03:01

Standard Lot#: 11514100A2

Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS

#53545

TUCSON PD

Operator: W. LOOMIS

#53545

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M

Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:03:54
Diagnostic Test	Pass	10:04:29
Air Blank	0.000	10:04:56
0.100 Cal Check	0.102	10:05:17
Air Blank	0.000	10:05:46
Subject Test	0.110*	10:08:54
Air Blank	0.000	10:09:26
Air Blank	0.000	10:09:53
Subject Test	XXX**	10:10:10
Air Blank	0.000	10:10:44

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

4 of 7  
wll  
Bh

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
02/04/2016 10:11:15

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:11:44
0.100 Cal Check	0.100	10:12:05
Air Blank	0.000	10:12:35

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
B4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL  
104 2.11.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 2/4/2016 TIME 09:54  
Intoxilyzer Serial # 80-002262 Location ODS

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE William Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
02/04/2016 09:54:25

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:54:55
0.100 Cal Check	0.101	09:55:15
Air Blank	0.000	09:55:45

Preventative Maintenance  
Performed -- Timer Reset

7 of 7  
WU  
B4