



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
3.23.2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/3/2016 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 2/4/2016 and 3/3/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


3/7/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 7th, day of March, 2016.

Elma Escarcega

Notary Public

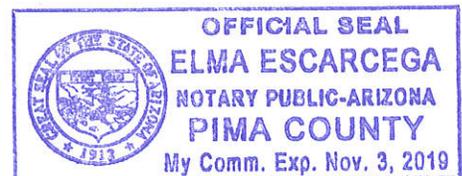

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
TS4 3-23-2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 3/3/2016 TIME 10:01

INTOXILYZER SERIAL # 80-002262 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Willie Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/03/2016

10:03:09

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:04:08 |
| Diagnostic Test | Pass | 10:04:43 |
| Air Blank | 0.000 | 10:05:11 |
| 0.100 Cal Check | 0.101 | 10:05:32 |
| Air Blank | 0.000 | 10:06:01 |
| Subject Test | 0.000 | 10:06:19 |
| Air Blank | 0.000 | 10:06:46 |
| Wait | ABT* | |
| Air Blank | 0.000 | 10:07:20 |

*Sequence Aborted

2 of 7
WLL
TSG

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/03/2016

10:07:26

3 of 7
WVW
TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:08:23 |
| Diagnostic Test | Pass | 10:08:58 |
| Air Blank | 0.000 | 10:09:25 |
| 0.100 Cal Check | RFI* | 10:09:39 |
| Air Blank | IPS** | 10:09:48 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/03/2016

10:09:56

4 of 7
WLL
T26

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:10:46 |
| Diagnostic Test | Pass | 10:11:21 |
| Air Blank | 0.000 | 10:11:49 |
| 0.100 Cal Check | 0.101 | 10:12:10 |
| Air Blank | 0.000 | 10:12:39 |
| Subject Test | 0.394* | 10:15:47 |
| Air Blank | 0.000 | 10:16:20 |
| Air Blank | 0.000 | 10:16:48 |
| Subject Test | XXX** | 10:17:09 |
| Air Blank | 0.000 | 10:17:43 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/03/2016 10:18:15

5 of 7
WLL
T34

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:18:45 |
| 0.100 Cal Check | 0.099 | 10:19:06 |
| Air Blank | 0.000 | 10:19:36 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
T54
3-23-2016
WLL

QA Specialist W. LOOMIS 53945 AGENCY Tucson Police Department
DATE 3/3/2016 TIME 10:01
Intoxilyzer Serial # 80-002262 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE William Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002262
Core Version: 8105.48
03/03/2016 10:01:36

7 of 7
WLL
TSG

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:02:05 |
| 0.100 Cal Check | 0.101 | 10:02:26 |
| Air Blank | 0.000 | 10:02:55 |

Preventative Maintenance
Performed -- Timer Reset