



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

Handwritten notes: T54, 6-6-2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/2016 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 4/28/2016 and 5/26/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on
The instrument was returned to service on

Handwritten signature: William L. Loomis

William L. Loomis


Handwritten date: 6/1/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 1 day of June, 2016.

Handwritten signature: A.M. Hernandez

Notary Public

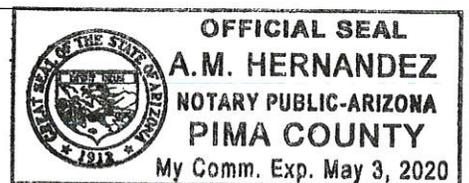

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
B4 6.6.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 5/26/2016 TIME 10:18

INTOXILYZER SERIAL # 80-002262 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Willie Loomis

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
05/26/2016

10:20:59

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:21:57
Diagnostic Test	Pass	10:22:32
Air Blank	0.000	10:22:59
0.100 Cal Check	0.101	10:23:20
Air Blank	0.000	10:23:50
Subject Test	0.000	10:24:08
Air Blank	0.000	10:24:35
Wait	ABT*	
Air Blank	0.000	10:25:09

\*Sequence Aborted

2 of 7  
WLL  
TS4

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
05/26/2016

10:25:14

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:26:12
Diagnostic Test	Pass	10:26:47
Air Blank	0.000	10:27:14
0.100 Cal Check	RFI*	10:27:28
Air Blank	IPS**	10:27:35

\*RFI Detect  
\*\*Improper Sample

3 of 7  
6 WLL  
TSL

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
05/26/2016

10:27:44

Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:28:39
Diagnostic Test	Pass	10:29:14
Air Blank	0.000	10:29:41
0.100 Cal Check	0.101	10:30:02
Air Blank	0.000	10:30:31
Subject Test	0.149*	10:33:39
Air Blank	0.000	10:34:11
Air Blank	0.000	10:34:39
Subject Test	XXX**	10:35:03
Air Blank	0.000	10:35:36

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

4 of 7  
WLL  
TSh

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
05/26/2016 10:35:57

-----  
Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:36:26
0.100 Cal Check	0.099	10:36:47
Air Blank	0.000	10:37:17

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
WLL  
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WB  
6.6.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 5/26/2016 TIME 10:18  
Intoxilyzer Serial # 80-002262 Location ODS

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002262  
Core Version: 8105.48  
05/26/2016 10:18:36

-----  
Standard Lot#: 06616100A1  
Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:19:10
0.100 Cal Check	0.100	10:19:31
Air Blank	0.000	10:20:00

Preventative Maintenance  
Performed -- Timer Reset

7 of 7  
0 WLL  
TS4