



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
6.12.07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-7-07 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 5-10-07 and 6-7-07
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

6-8-2007

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8, day of
June, 2007.

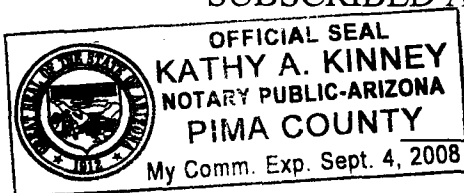

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley 34517 AGENCY Tucson Police DepartmentDATE 6-7-07 TIME 1104INTOXILYZER SERIAL # 80-002264 LOCATION Traffic

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 ACInstrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Clare Conley

1 of 9 cu cu
1 of 8 cu
B4612-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
06/07/2007

11:53:53

Standard Lot#: 629202I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:54:49
Diagnostic Test	Pass	11:55:17
Air Blank	0.000	11:55:44
0.100 Cal Check	0.101	11:56:03
Air Blank	0.000	11:56:32
Subject Test	0.000	11:56:49
Air Blank	ABT*	11:56:59
Air Blank	0.000	11:57:27

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 8 cu
BY
6-12-07

3 of 8 cu

By
6-12-07

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

06/07/2007

12:00:45

Standard Lot#: 629202I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:03:06
Diagnostic Test	Pass	12:03:34
Air Blank	0.000	12:04:01
0.100 Cal Check	RFI*	12:04:13
Air Blank	0.000	12:04:44

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 Cn

TS4
6-12-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
06/07/2007

11:57:34

Standard Lot#: 629202I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:58:29
Diagnostic Test	Pass	11:58:57
Air Blank	0.000	11:59:24
0.100 Cal Check	0.008*	11:59:43
Air Blank	IPS**	12:00:04

*Calibration Check Out of Tolerance

**Improper Sample

*Blew during calibration check and then
blew during subsequent air blank.*

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 8 Cu

B4
6-12-07

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

06/07/2007

12:04:51

Standard Lot#: 629202I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:05:44
Diagnostic Test	Pass	12:06:12
Air Blank	0.000	12:06:39
0.100 Cal Check	0.102	12:06:58
Air Blank	0.000	12:07:27
Subject Test	XXX*	12:10:24
Air Blank	0.000	12:10:59

*Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

06/07/2007

12:11:39

Standard Lot#: 629202I

Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:12:07
Cal Check	0.101	12:12:25
Air Blank	0.000	12:12:55

Clare Conley

Operator's Signature

6 of 8 cu

B4
6.12.07

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley 34517 AGENCY Tucson Police Department
DATE 6-7-07 TIME 1104
Intoxilyzer Serial # 80-00 2264 Location Traffic

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

06/07/2007

11:52:27

Standard Lot#: 629202I

Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:52:55
Cal Check	0.102	11:53:14
Air Blank	0.000	11:53:43

Clare Conley

Operator's Signature

CM
~~8/8 CM~~
8/8 CM
By
6.11.07