



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

Ch
2-20-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/14/08 Cal 1 0.103 AC Cal 2 0.102 AC

Between the dates of 1/16/08 and 2/14/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

2/15/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15, day of February, 2008.

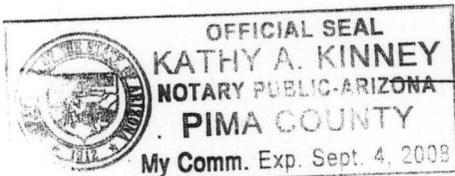

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

Ch

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

2-20-08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/14/08 TIME 10:51

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES NO

COMMENTS none

SIGNATURE Sarah Cione

2 of 7 sc

CU
2-20-08

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-102264
02/14/2008 10:54:51

Standard Lot#: 7254021
Last Changed By: J. HAWKINS #42401

DPS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SCAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:55:59
Diagnostic Test	Pass	10:56:27
Air Blank	0.000	10:56:54
0.100 Cal Check	0.103	10:57:13
Air Blank	0.000	10:57:42
Subject Test	0.000	10:58:00
Air Blank	0.000	10:58:28
Five Minute wait	ABT*	
Air Blank	0.000	10:59:46

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7 sc

CU

2-20-08

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
02/14/2008 11:09:01

Standard Lot#: 7254021
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:09:42
Air Blank	IPS**	11:10:03

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7sc

du

2-20-08

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002254
02/14/2006 11:00:31

Standard Lock: 7254021
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:01:44
Diagnostic Test	Pass	11:02:13
Air Blank	0.000	11:02:40
0.100 Cal Check	0.103	11:02:59
Air Blank	0.000	11:03:29
Subject Test	0.500*	11:06:37
Air Blank	0.000	11:07:10
Air Blank	0.000	11:07:37
Subject Test	XXX**	11:08:06
Air Blank	0.000	11:08:39

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7sc
CU
2-20-08

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 60-002264
02/14/2008 11:10:26

Standard Lot#: 7254021
Last Changed By: J. HAWKINS #4240.

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:10:53
Cal Check	0.102	11:11:12
Air Blank	0.000	11:11:41

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7sc
U
2-20-08

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 2/14/08 TIME 10:51

Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

02/14/2008

10:51:49

7 of 7sc
cu
2-22-08

Standard Lot#: 7254021

Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:52:16
Cal Check	0.103	10:52:35
Air Blank	0.000	10:53:04

Sarah Cione

Operator's Signature