



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-802264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

BY  
6-6-2008

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-5-08 Cal 1 0.102 AC Cal 2 0.102 AC

Between the dates of 5-8-08 and 6-5-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A  
The instrument was returned to service on N/A


Clare Conley

Clare L. Conley

6-6-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 6, day of June, 2008.

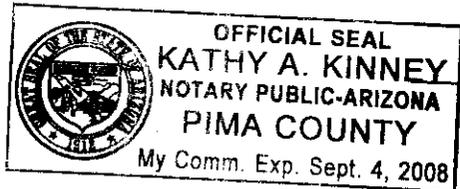

Kathy A Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1567  
Cu  
TSY  
6-6-2008

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 6-5-08 TIME 0956

INTOXILYZER SERIAL # 80 - 002264 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 102 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
06/05/2008

09:57:55

2 of 7  
cu  
PH

Standard Lot#: 735401G  
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:59:04
Diagnostic Test	Pass	09:59:32
Air Blank	0.000	09:59:59
0.100 Cal Check	0.101	10:00:18
Air Blank	0.000	10:00:47
Subject Test	0.000	10:01:06
Air Blank	ABT*	10:01:13
Air Blank	0.000	10:01:41

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
06/05/2008

10:01:51

Standard Lot#: 735401G  
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:02:44
Air Blank	IPS**	10:03:01

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

3587  
CW

154

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
06/05/2008

10:03:12

4687  
cu  
SM

Standard Lot#: 735401G  
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:04:10
Diagnostic Test	Pass	10:04:38
Air Blank	0.000	10:05:05
0.100 Cal Check	0.102	10:05:24
Air Blank	0.000	10:05:53
Subject Test	0.079*	10:09:02
Air Blank	0.000	10:09:33
Air Blank	0.000	10:10:00
Subject Test	XXX**	10:10:24
Air Blank	0.000	10:10:58

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002264  
06/05/2008

10:12:33

-----  
Standard Lot#: 735401G  
Last Changed By: J. HAWKINS #42401  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:13:01
Cal Check	0.102	10:13:19
Air Blank	0.000	10:13:49

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

5067  
cu  
154

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

647  
CM  
TSC  
6-6-08

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 6-5-08 TIME 0956  
Intoxilyzer Serial # 80-002264 Location DUI HQ ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

06/05/2008

09:56:20

7067  
cu

-----  
Standard Lot#: 735401G

Last Changed By: J. HAWKINS #42401  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

TS4

Test	g/210L	Time
Air Blank	0.000	09:56:47
Cal Check	0.102	09:57:06
Air Blank	0.000	09:57:36

*Clare Conley*

Operator's Signature