



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

U
7-10-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/3/08 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 6/5/08 and 7/3/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione


Sarah L. Cione

7/3/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3, day of July, 2008.

Kathy A Kinney
Notary Public

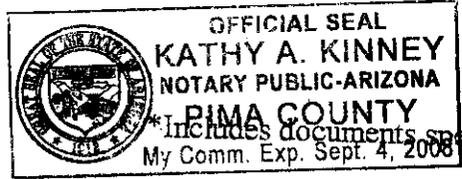

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

DU
7-10-08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/3/08 TIME 9:24

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7 sc

cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
07/03/2008 09:26:13

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:26:57
Diagnostic Test	Pass	09:27:24
Air Blank	0.000	09:27:51
0.100 Cal Check	0.101	09:28:10
Air Blank	0.000	09:28:39
Subject Test	0.000	09:29:06
Air Blank	0.000	09:29:34
Five Minute Wait	ABT*	
Air Blank	0.000	09:30:10

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
07/03/2008

09:30:17

3 of 75c

cu

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:30:55
Air Blank	IPS**	09:31:05

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
07/03/2008

09:31:39

4 of 75c

cu

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:32:34
Diagnostic Test	Pass	09:33:01
Air Blank	0.000	09:33:29
0.100 Cal Check	0.100	09:33:48
Air Blank	0.000	09:34:17
Subject Test	0.073*	09:37:25
Air Blank	0.000	09:37:57
Air Blank	0.000	09:38:25
Subject Test	XXX**	09:38:48
Air Blank	0.000	09:39:21

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
07/03/2008 09:39:38

5 of 7 sc

Ch

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Test g/210L Time

Air Blank 0.000 09:40:05
Cal Check 0.100 09:40:24
Air Blank 0.000 09:40:54



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

CU
7-10-08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Ciene AGENCY Tucson Police Department
DATE 7/3/08 TIME 9:24
Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Ciene

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
07/03/2008 09:24:14

7 of 7sc

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

ck

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:24:42
Cal Check	0.101	09:25:00
Air Blank	0.000	09:25:30

Sarah Cione
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2 sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

OK
7-10-08

Date: 6/26/08
QAS: HAWKINS 42401
Intoxilyzer Serial Number: 80002264
New Lot Number: 735401G
New Expiration Date: 28 Dec 2009

(S)
(S)
(S)

1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
2. Remove the old gas tank and install the new one.
3. Hit the "ESC" key twice and swipe the QAS card.
4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:

Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."

(S)
(S)
(S)
(S)
()

5. Hit the "ESC" key twice.
6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
7. Air blank completed.
8. Calibration check completed. Test results 0. 101 AC.
9. Air blank completed.
10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

[Handwritten Signature]

2 of 25c

CM

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
06/26/2008 22:59:50

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

OPS: J. HAWKINS #42401
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	23:00:17
Cal Check	0.101	23:00:36
Air Blank	0.000	23:01:05


Operator's Signature