



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

By
8-29-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-28-08 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 7-31-08 and 8-28-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A

The instrument was returned to service on N/A


Clare Conley

Clare L. Conley

8-29-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 29th, day of August, 2008.

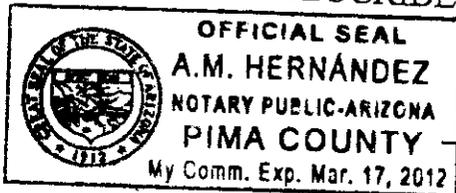

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1067
cu
By
8-29-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 8-28-08 TIME 1124

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ Team 2

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/28/2008

11:25:50

287
cu

104

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:26:37
Diagnostic Test	Pass	11:27:05
Air Blank	0.000	11:27:32
0.100 Cal Check	0.100	11:27:51
Air Blank	0.000	11:28:20
Subject Test	0.000	11:28:37
Air Blank	ABT*	11:28:43
Air Blank	0.000	11:29:12

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/28/2008

11:29:27

387
cu

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:30:02
Air Blank	IPS**	11:30:11

*RFI Detect

**Improper Sample

BH

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/28/2008

11:31:01

487
cu

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

134

Test	g/210L	Time
Air Blank	0.000	11:31:57
Diagnostic Test	Pass	11:32:24
Air Blank	0.000	11:32:51
0.100 Cal Check	0.100	11:33:10
Air Blank	0.000	11:33:40
Subject Test	0.132*	11:36:48
Air Blank	0.000	11:37:20
Air Blank	0.000	11:37:48
Subject Test	XXX**	11:38:07
Air Blank	0.000	11:38:41

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

754
8-29-08
657
cu

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 8-28-08 TIME 1124
Intoxilyzer Serial # 80-002264 Location DUI HQ Team 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
08/28/2008 11:24:29

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:24:56
Cal Check	0.100	11:25:15
Air Blank	0.000	11:25:44



Operator's Signature

7507
Cu
7504