



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

SC
11/21/08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11-20-08 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 10-23-08 and 11-20-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


11-20-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 20th, day of November, 2008.

A.A. Hernandez

Notary Public

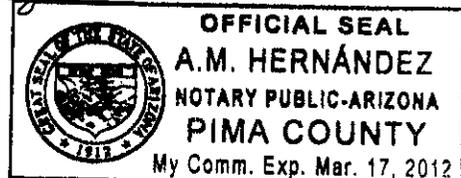

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1577 cu

SC
11/21/08

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 11-20-08 TIME 0930

INTOXILYZER SERIAL # 80-002264 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
11/20/2008

09:32:19

2587
lu
sc

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:33:07
Diagnostic Test	Pass	09:33:34
Air Blank	0.000	09:34:01
0.100 Cal Check	0.103	09:34:20
Air Blank	0.000	09:34:49
Subject Test	0.000	09:35:09
Air Blank	ABT*	09:35:15
Air Blank	0.000	09:35:43

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
11/20/2008

09:36:40

3567
LW
SC

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:37:09
Air Blank	IPS**	09:37:17

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
11/20/2008

09:38:10

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:39:12
Diagnostic Test	Pass	09:39:39
Air Blank	0.000	09:40:07
0.100 Cal Check	0.102	09:40:25
Air Blank	0.000	09:40:55
Subject Test	0.167*	09:44:03
Air Blank	0.000	09:44:35
Air Blank	0.000	09:45:03
Subject Test	XXX**	09:45:29
Air Blank	0.000	09:46:03

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

487
lw
SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
11/20/2008

09:46:22

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:46:49
Cal Check	0.101	09:47:08
Air Blank	0.000	09:47:37

Clare Conley

Operator's Signature

5067
Cu
SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6587
CW
SC

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 11-20-08 TIME 0930
Intoxilyzer Serial # 80-002264 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
11/20/2008 09:30:52

7567
CW
SC

Standard Lot#: 735401G
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:31:20
Cal Check	0.102	09:31:39
Air Blank	0.000	09:32:08

Clare Conley

Operator's Signature