



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
 3-17-09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-12-09 Cal 1 0.102 AC Cal 2 0.099 AC

Between the dates of 2-12-09 and 3-12-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
 The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


3-17-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 17, day of March, 2009.

Amy Hernandez

Notary Public

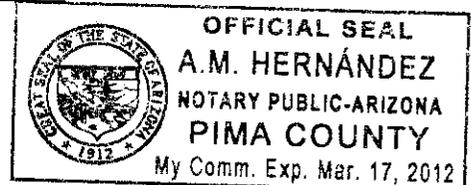

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1587

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

cu

STANDARD QUALITY ASSURANCE PROCEDURE

154
3/17/09

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 3-12-09 TIME 0906

INTOXILYZER SERIAL # 80-002264 LOCATION ODW DWI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
03/12/2009

09:10:34

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:11:32
Diagnostic Test	Pass	09:12:00
Air Blank	0.000	09:12:27
0.100 Cal Check	0.101	09:12:46
Air Blank	0.000	09:13:15
Subject Test	0.000	09:13:57
Air Blank	ABT*	09:14:04
Air Blank	0.000	09:14:32

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

2577
lu
154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
03/12/2009 09:15:02

3567
en
154

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:15:56
Air Blank	IPS**	09:16:06

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
03/12/2009

09:16:13

487
en
154

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:16:58
Diagnostic Test	Pass	09:17:25
Air Blank	0.000	09:17:53
0.100 Cal Check	0.101	09:18:11
Air Blank	0.000	09:18:41
Subject Test	0.499*	09:21:49
Air Blank	0.000	09:22:24
Air Blank	0.000	09:22:52
Subject Test	XXX**	09:23:32
Air Blank	0.000	09:24:06

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/12/2009 09:24:31

5-67
cu
TS4

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:24:58
Cal Check	0.099	09:25:17
Air Blank	0.000	09:25:47

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6887
lu

TSK
3-17-09

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 3-12-09 TIME 0906
Intoxilyzer Serial # 80-002264 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/12/2009 09:06:51

7-67
lu

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

154

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:07:18
Cal Check	0.102	09:07:37
Air Blank	0.000	09:08:06

Clare Conley

Operator's Signature