



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-00 2264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

Handwritten note: T54 4.10.09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-9-09 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 3-12-09 and 4-9-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A. The instrument was returned to service on N/A.

Signature of Clare Conley

Clare L. Conley


4-9-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 9th day of April, 2009.

Signature of Notary Public

Notary Public

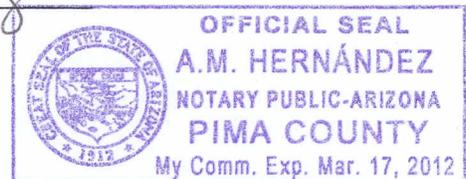

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1577
an

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSH
4.10.09

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 4-9-09 TIME 0752

INTOXILYZER SERIAL # 80-002264 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
04/09/2009

07:53:57

287
On
TS4

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:54:43
Diagnostic Test	Pass	07:55:11
Air Blank	0.000	07:55:38
0.100 Cal Check	0.101	07:55:57
Air Blank	0.000	07:56:26
Subject Test	0.000	07:57:44
Air Blank	ABT*	07:57:50
Air Blank	0.000	07:58:19

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
04/09/2009 07:58:23

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:58:53
Air Blank	IPS**	07:59:01

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

3587
cu
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
04/09/2009

07:59:33

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:00:21
Diagnostic Test	Pass	08:00:48
Air Blank	0.000	08:01:16
0.100 Cal Check	0.101	08:01:34
Air Blank	0.000	08:02:04
Subject Test	0.106*	08:05:12
Air Blank	0.000	08:05:47
Air Blank	0.000	08:06:15
Subject Test	XXX**	08:06:53
Air Blank	0.000	08:07:27

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

4587
Or
TS4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

04/09/2009

08:07:45

Standard Lot#: 811201I

Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:08:12
Cal Check	0.101	08:08:31
Air Blank	0.000	08:09:00

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

b87
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TSY
4.10.09

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 4-9-09 TIME 0752
Intoxilyzer Serial # 80-002264 Location ODW DUI HA

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
04/09/2009 07:52:29

787
CM

Standard Lot#: 811201I
Last Changed By: J. HAWKINS #42401

154

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:52:56
Cal Check	0.102	07:53:15
Air Blank	0.000	07:53:44

C. Conley

Operator's Signature