



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
6-26-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/25/09 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 6/4/09 and 6/25/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/25/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 25<sup>th</sup> day of June, 20 09.

A.M. Hernandez

Notary Public

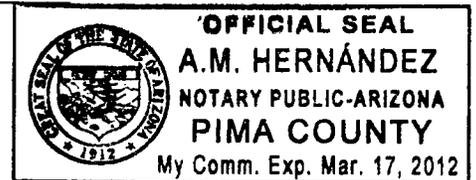

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TSY  
626-09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 6/25/09 TIME 9:22  
INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7sc

BY

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
06/25/2009 09:23:55

Standard Lot#: 730901I  
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:24:40
Diagnostic Test	Pass	09:25:07
Air Blank	0.000	09:25:35
0.100 Cal Check	0.100	09:25:54
Air Blank	0.000	09:26:23
Subject Test	0.000	09:26:42
Air Blank	0.000	09:27:09
Five Minute Wait	ABT*	
Air Blank	0.000	09:27:47

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7sc

ASU

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
06/25/2009 09:29:22

---

Standard Lot#: 730901I  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	09:29:51
Air Blank	IPS**	09:30:01

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>

154

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
06/25/2009 09:30:08

Standard Lot#: 730901I  
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:31:10
Diagnostic Test	Pass	09:31:37
Air Blank	0.000	09:32:05
0.100 Cal Check	0.101	09:32:23
Air Blank	0.000	09:32:53
Subject Test	0.122*	09:36:01
Air Blank	0.000	09:36:33
Air Blank	0.000	09:37:00
Subject Test	XXX**	09:37:19
Air Blank	0.000	09:37:52

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc

TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002264  
06/25/2009 09:38:08

-----  
Standard Lot#: 730901I  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:38:36
Cal Check	0.100	09:38:54
Air Blank	0.000	09:39:24

*Sarah Cione*

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSY  
6.26.09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/25/09 TIME 9:22  
Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7sc

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002264  
06/25/2009 09:22:12

By

-----  
Standard Lot#: 730901I  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:22:40
Cal Check	0.101	09:22:59
Air Blank	0.000	09:23:28

Sarah Cione  
Operator's Signature