



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
8.31.09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/20/09 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 07/23/09 and 08/20/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/31/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 2009.

A.M. Hernandez

Notary Public

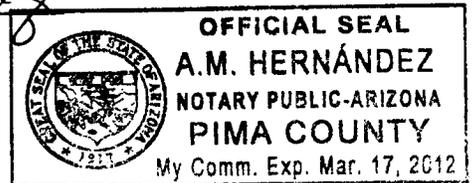

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-1041A

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

By
8.31.2009

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department
DATE 08/20/09 TIME 11:40
INTOXILYZER SERIAL # 80-002264 LOCATION DUT SQUAD HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Cheryl Straub

287
C8

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/20/2009

11:43:27

B4

Standard Lot#: 730901I
Last Changed By: J. HAWKINS #42401

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:44:55
Diagnostic Test	Pass	11:45:22
Air Blank	0.000	11:45:50
0.100 Cal Check	0.100	11:46:08
Air Blank	0.000	11:46:38
Subject Test	0.000	11:47:02
Air Blank	ABT*	11:47:11
Air Blank	0.000	11:47:39

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

387
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/20/2009 11:48:02

BY

Standard Lot#: 730901I
Last Changed By: J. HAWKINS #42401

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	α/210L	Time
Air Blank	RFI*	11:49:10
Air Blank	IPS**	11:49:22

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

489
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/20/2009

11:49:45

754

Standard Lot#: 730901I
Last Changed By: J. HAWKINS #42401

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:50:40
Diagnostic Test	Pass	11:51:07
Air Blank	0.000	11:51:35
0.100 Cal Check	0.100	11:51:54
Air Blank	0.000	11:52:23
Subject Test	0.179*	11:55:31
Air Blank	0.000	11:56:04
Air Blank	0.000	11:56:32
Subject Test	XXX**	11:56:46
Air Blank	0.000	11:57:20

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

597
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

08/20/2009

11:58:12

PSY


Standard Lot#: 730901I

Last Changed By: J. HAWKINS #42401

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:58:40
Cal Check	0.099	11:58:58
Air Blank	0.000	11:59:28

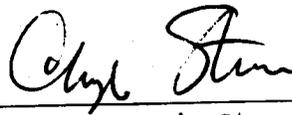

Operator's Signature

687 CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

By 8.31.09

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 08/20/09 TIME 11:40
Intoxilyzer Serial # 80-002264 Location DUI SQUAD HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787
es

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
08/20/2009 11:40:49

B4

Standard Lot#: 730901I
Last Changed By: J. HAWKINS #42401

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:41:17
Cal Check	0.100	11:41:35
Air Blank	0.000	11:42:05



Operator's Signature

172
cs



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 8/23/09
QAS: HAWKINS 42401
Intoxilyzer Serial Number: 80-002264
New Lot Number: 826001I
New Expiration Date: 25 Sep 2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE [Signature]

117
CS

INTOKI 228-3111
Location: TUCSON POLICE DEPT
Serial Number: 80-102264
08/23/2009 21:25:38

Standard Lot#: 8260011
Last Changed By: J. HAWKINS #42401

DAS: J. HAWKINS #42401
TUCSON PD

Test	q/210L	Time
Ans. Blank	0.000	21:26:05
Cal. Check	0.000	21:26:24
Ans. Blank	0.000	21:26:53

Operator: J. Hawkins