



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
5-5-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/29/10 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 4/1/10 and 4/29/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/30/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30th, day of April, 2010.

A M Hernandez

Notary Public

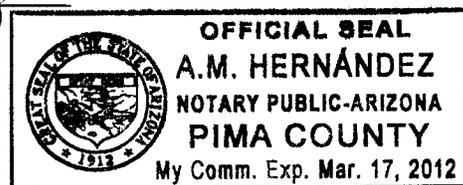

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TJH
5-5-2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/29/10 TIME 959

INTOXILYZER SERIAL # 80-002264 LOCATION DUT HQ - C12

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
04/29/2010 10:01:40

2 of 7 sc

TS4

Standard Lot#: 642845
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:02:23
Diagnostic Test	Pass	10:02:51
Air Blank	0.000	10:03:18
0.100 Cal Check	0.101	10:03:37
Air Blank	0.000	10:04:06
Subject Test	0.000	10:05:05
Air Blank	0.000	10:05:33
Five Minute Wait	ABT*	
Air Blank	0.000	10:07:11

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
04/29/2010 10:07:19

3 of 7 sc
BH

Standard Lot#: 642845
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:08:03
Diagnostic Test	Pass	10:08:31
Air Blank	RFI*	10:08:39
Air Blank	IPS**	10:08:47

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}
ADY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
04/29/2010 10:08:53

Standard Lot#: 642845
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:09:39
Diagnostic Test	Pass	10:10:06
Air Blank	0.000	10:10:33
0.100 Cal Check	0.100	10:10:52
Air Blank	0.000	10:11:22
Subject Test	0.047*	10:14:30
Air Blank	0.000	10:15:01
Air Blank	0.000	10:15:28
Subject Test	XXX**	10:15:46
Air Blank	0.000	10:16:19

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 9000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
04/29/2010 10:16:36

5 of 7 SC
104

Standard Lot#: 642845
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:17:03
Cal Check	0.101	10:17:22
Air Blank	0.000	10:17:52

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
5-5-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 4/29/10 TIME 959
Intoxilyzer Serial # 80-002264 Location C12-DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
04/29/2010 09:59:46

7 of 7 sc

124

Standard Lot#: 642845
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:00:13
Cal Check	0.101	10:00:32
Air Blank	0.000	10:01:01

Sarah Cione
Operator's Signature