



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
6.7-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/10 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 4/29/10 and 5/26/10
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

6/3/10


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3rd, day of June, 2010.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

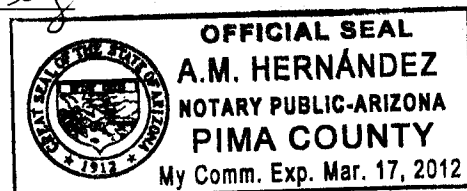

1 of 8 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B4
6-7-2010QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 5/26/10 TIME 852INTOXILYZER SERIAL # 80-002264 LOCATION Team 2 - DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Sarah Cione

2 of 8_{sc}

ph

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
05/26/2010

08:54:24

Standard Lot#: 642845
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:55:50
Diagnostic Test	Pass	08:56:18
Air Blank	0.000	08:56:45
0.100 Cal Check	0.101	08:57:04
Air Blank	0.000	08:57:34
Subject Test	0.000	08:58:27
Air Blank	0.000	08:58:54
Five Minute Wait	ABT*	
Air Blank	0.000	08:59:44

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8_{sc}

B4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

05/26/2010

08:59:49

Standard Lot#: 642845

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:00:20
Air Blank	IPS**	09:00:28

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

05/26/2010

09:00:59

4 of 8_{sc}

BH

Standard Lot#: 642845

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:01:57
Diagnostic Test	Pass	09:02:24
Air Blank	0.000	09:02:52
0.100 Cal Check	0.102	09:03:11
Air Blank	0.000	09:03:40
Subject Test	RFI*	09:05:25
Air Blank	0.000	09:05:57

*RFI Detect

Not a Successfully Completed Test Sequence

5 of 8_{sc}

189

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

05/26/2010

09:06:03

Standard Lot#: 642845

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:07:12
Diagnostic Test	Pass	09:07:39
Air Blank	0.000	09:08:07
0.100 Cal Check	0.101	09:08:26
Air Blank	0.000	09:08:55
Subject Test	0.024*	09:12:03
Air Blank	0.000	09:12:34
Air Blank	0.000	09:13:01
Subject Test	XXX**	09:13:53
Air Blank	0.000	09:14:25

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

6 of 8 sc

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

05/26/2010

09:17:47

BH

Standard Lot#: 642845

Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:18:14
Cal Check	0.101	09:18:33
Air Blank	0.000	09:19:03

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8_{sc}

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000BY
6.7.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/26/10 TIME 852
Intoxilyzer Serial # 80-002264 Location Team 2 - DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

05/26/2010

08:52:32

Standard Lot#: 642845

Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:52:59
Cal Check	0.101	08:53:18
Air Blank	0.000	08:53:48

Sarah Cione

Operator's Signature

8 of 8 sc

B4