



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

B4
6-25-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/24/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 5/26/10 and 6/24/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

6/24/10


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 24, day of June, 2010.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

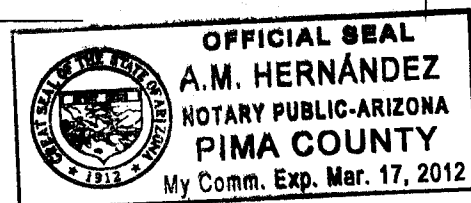

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY
6.25.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/24/10 TIME 1332

INTOXILYZER SERIAL # 80-002264 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0.100 AC

Result: 0.100 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Sarah Cione

2 of 7_{sc}

184

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

06/24/2010

13:32:25

Standard Lot#: 642845

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:33:09
Diagnostic Test	Pass	13:33:36
Air Blank	0.000	13:34:04
0.100 Cal Check	0.101	13:34:23
Air Blank	0.000	13:34:52
Subject Test	0.000	13:35:37
Air Blank	0.000	13:36:04
Five Minute Wait	ABT*	
Air Blank	0.000	13:37:04

*Sequence Aborted

Not a Successfully Completed Test Sequence

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
06/24/2010

13:37:34

Standard Lot#: 642845
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:38:37
Diagnostic Test	Pass	13:39:05
Air Blank	0.000	13:39:32
0.100 Cal Check	0.101	13:39:51
Air Blank	0.000	13:40:20
Subject Test	RFI*	13:40:48
Air Blank	IPS**	13:40:57

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

06/24/2010

13:41:14

BM

Standard Lot#: 642845

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:41:59
Diagnostic Test	Pass	13:42:26
Air Blank	0.000	13:42:54
0.100 Cal Check	0.101	13:43:12
Air Blank	0.000	13:43:42
Subject Test	0.114*	13:46:50
Air Blank	0.000	13:47:23
Air Blank	0.000	13:47:50
Subject Test	XXX**	13:48:24
Air Blank	0.000	13:48:57

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7_x

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

06/24/2010

13:49:12

TJH


Standard Lot#: 642845

Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:49:40
Cal Check	0.100	13:49:58
Air Blank	0.000	13:50:28

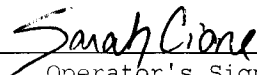

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/24/10 TIME 1332
Intoxilyzer Serial # 80-002264 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

06/24/2010

13:30:51

Standard Lot#: 642845

Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:31:18
Cal Check	0.101	13:31:37
Air Blank	0.000	13:32:07

Sarah Cione
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

1 of 1 sc

T34
7-6-2010

Date: 6/24/10

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002264

Reason For Service: Officer Hawkins left note that there was a temp. regulation failure. When I turned intox on, breath hose would not heat up. I brought intox to the lab to troubleshoot.

Date Out of Service: N/A

Calibration Check completed prior to removal from service

Date: -

SQAP completed prior to removal from service

Date: -

Corrective Action:

☐ Sent to CMI, Inc.

Removed breath hose + replaced with a new breath hose. Was able to complete SQAP + did not have to remove instrument from service.