



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

TS9
 8-23-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/19/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 7/22/10 and 8/19/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
 The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


8/23/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 23, day of August, 2010.

A.M. Hernandez

Notary Public

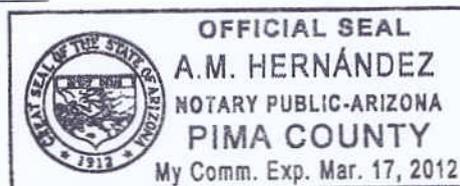

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
8-23-2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 8/19/10 TIME 915

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

BS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/19/2010 09:20:13

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:20:59
Diagnostic Test	Pass	09:21:26
Air Blank	0.000	09:21:53
0.100 Cal Check	0.100	09:22:12
Air Blank	0.000	09:22:42
Subject Test	0.000	09:23:38
Air Blank	0.000	09:24:05
Five Minute Wait	ABT*	
Air Blank	0.000	09:24:57

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc
154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/19/2010 09:25:32

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:26:07
Air Blank	IPS**	09:26:15

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc
154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
08/19/2010

09:27:57

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:28:40
Diagnostic Test	Pass	09:29:08
Air Blank	0.000	09:29:35
0.100 Cal Check	0.100	09:29:54
Air Blank	0.000	09:30:24
Subject Test	0.189*	09:33:32
Air Blank	0.000	09:34:04
Air Blank	0.000	09:34:31
Subject Test	XXX**	09:35:09
Air Blank	0.000	09:35:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
08/19/2010 09:36:12

5 of 7 SC
BS4

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:36:39
Cal Check	0.100	09:36:58
Air Blank	0.000	09:37:27

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
8.23.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 8/19/10 TIME 915
Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
08/19/2010 09:15:26

7 of 7 SC
154

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:15:53
Cal Check	0.101	09:16:12
Air Blank	0.000	09:16:42


Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2 sc

TS4
8-23-2010

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 07/26/10
 QAS: HAWKINS 424/01
 Intoxilyzer Serial Number: 80-002264
 New Lot Number: 719919
 New Expiration Date: 02 February 2012

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE _____


2 of 2 sc
TX

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
07/26/2010 23:44:10

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

DAS: J. HAWKINS #42401
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	23:44:37
Cal Check	0.100	23:44:56
Air Blank	0.000	23:45:26

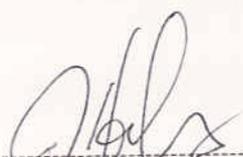
Operator's Signature