



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TSY
10-18-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 9/16/10 and 10/14/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


10/14/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14, day of October, 2010.

A.M. Hernandez

Notary Public

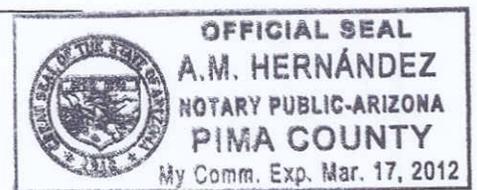

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8_{sc}

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

By
10-18-2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/10 TIME 838

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS Interrupted while performing operational tests - performed alcohol free
subject test and sequence abort twice

SIGNATURE Sarah Cione

2 of 8_{sc}

B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
10/14/2010 08:39:51

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:40:35
Diagnostic Test	Pass	08:41:03
Air Blank	0.000	08:41:30
0.100 Cal Check	0.100	08:41:49
Air Blank	0.000	08:42:18
Subject Test	0.000	08:42:52
Air Blank	0.000	08:43:19
Five Minute Wait	ABT*	
Air Blank	0.000	08:44:04

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 8 sc

124

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
10/14/2010 08:44:37

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:45:10
Air Blank	IPS**	08:45:17

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
10/14/2010 09:07:20

4 of 8 SC

134

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:08:08
Diagnostic Test	Pass	09:08:36
Air Blank	0.000	09:09:03
0.100 Cal Check	0.100	09:09:22
Air Blank	0.000	09:09:51
Subject Test	0.000	09:10:40
Air Blank	0.000	09:11:08
Five Minute Wait	ABT*	
Air Blank	0.000	09:11:50

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 8 sc

BY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
10/14/2010 09:12:00

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:12:48
Diagnostic Test	Pass	09:13:15
Air Blank	0.000	09:13:43
0.100 Cal Check	0.100	09:14:01
Air Blank	0.000	09:14:31
Subject Test	0.058*	09:17:39
Air Blank	0.000	09:18:11
Air Blank	0.000	09:18:38
Subject Test	XXX**	09:18:56
Air Blank	0.000	09:19:29

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

10/14/2010

09:19:59

6 of 8 sc

TS4

Standard Lot#: 719919

Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:20:26
Cal Check	0.100	09:20:45
Air Blank	0.000	09:21:15

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

By
10-18-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/14/10 TIME 838
Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

8 of 8 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
10/14/2010 08:38:13

Standard Lot#: 719919
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

154

Test	g/210L	Time
Air Blank	0.000	08:38:41
Cal Check	0.101	08:39:00
Air Blank	0.000	08:39:29

Sarah Cione
Operator's Signature