



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
3-21-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-3-11 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 2-3-11 and 3-3-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


3-16-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 16th, day of March, 20 11.

Isabel C. Trujillo

Notary Public

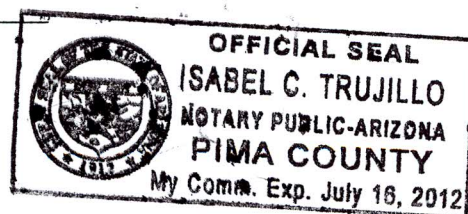

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1567
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSK
3.21.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 3-3-11 TIME 0743
INTOXILYZER SERIAL # 80-002264 LOCATION ODNDU HO

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS None.

SIGNATURE Clare Conley

2017
eu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/03/2011 07:45:29

TS4

Standard Lot#: 20910100A1
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:46:21
Diagnostic Test	Pass	07:46:49
Air Blank	0.000	07:47:16
0.100 Cal Check	0.101	07:47:35
Air Blank	0.000	07:48:04
Subject Test	0.000	07:48:35
Air Blank	ABT*	07:48:41
Air Blank	0.000	07:49:10

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3-67
cm

By

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/03/2011 07:50:55

Standard Lot#: 20910100A1
Last Changed By: J. HAWKINS #42401

OAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:51:24
Air Blank	IPS**	07:51:33

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/03/2011 07:52:16

Standard Lot#: 20910100A1
Last Changed By: J. HAWKINS #42401

OAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SOAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:53:03
Diagnostic Test	Pass	07:53:32
Air Blank	0.000	07:54:00
0.100 Cal Check	0.100	07:54:18
Air Blank	0.000	07:54:48
Subject Test	0.174*	07:57:56
Air Blank	0.000	07:58:28
Air Blank	0.000	07:58:55
Subject Test	XXX**	07:59:24
Air Blank	AMB***	07:59:58

*Deficient Sample

**Invalid Sample

Begin new deprivation period

***Ambient Fail

Not a Successfully
Completed
Test Sequence

5067
cu

ph

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/03/2011 08:01:11

Standard Lot#: 20910100A1
Last Changed By: J. HAWKINS #42401

OAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:01:39
Cal Check	0.099	08:01:57
Air Blank	0.000	08:02:27

Clare Conley
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

b-67
Cu
TS4
3.21.2011

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 3-3-11 TIME 0743
Intoxilyzer Serial # 80-002264 Location ODW Dai Hb

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

03/03/2011

07:43:57

Standard Lot#: 20910100A1

Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:44:25
Cal Check	0.101	07:44:43
Air Blank	0.000	07:45:13

Clare Conley

Operator's Signature

7567

CU

754