



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TSY
5-2-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-28-11 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 3-31-11 and 4-28-11
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


4-28-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 28, day of April, 20 11.

Isabel C. Trujillo

Notary Public

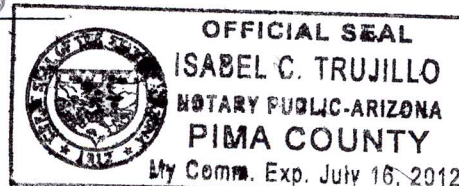

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1567
an

TS4
5-2-2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 4-28-11 TIME 0932

INTOXILYZER SERIAL # 80-002264 LOCATION DDW DM HQ

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(☒) 1. Clock time check.

(☒) 2. Date check.

OPERATIONAL TESTS

(☒) 1. Alcohol-free subject test result 0. 000 AC.

(☒) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(☒) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(☒) 4. Standard Calibration Check standard 0. 100 AC

Result: 0. 100 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS None.

SIGNATURE Clare Conley

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
04/28/2011

09:33:30

2567
CW
TD4

Standard Lot#: 20910100A1
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:34:20
Diagnostic Test	Pass	09:34:47
Air Blank	0.000	09:35:15
0.100 Cal Check	0.102	09:35:34
Air Blank	0.000	09:36:03
Subject Test	0.000	09:36:29
Air Blank	ABT*	09:36:36
Air Blank	0.000	09:37:04

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
04/28/2011

09:39:40

Standard Lot#: 20910100A1
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TES, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:40:20
Air Blank	IPS**	09:40:28

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

04/28/2011

09:41:05

Standard Lot#: 20910100A1

Last Changed By: J. HAWKINS

#42401

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:41:54
Diagnostic Test	Pass	09:42:22
Air Blank	0.000	09:42:49
0.100 Cal Check	0.101	09:43:08
Air Blank	0.000	09:43:37
Subject Test	0.211*	09:46:46
Air Blank	0.000	09:47:18
Air Blank	0.000	09:47:45
Subject Test	XXX**	09:47:59
Air Blank	0.000	09:48:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

04/28/2011

09:49:04

Standard Lot#: 20910100A1

Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:49:32
Cal Check	0.100	09:49:50
Air Blank	0.000	09:50:20



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 4-28-11 TIME 0932
Intoxilyzer Serial # 80-002264 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
04/28/2011 09:32:10

7587
m

Standard Lot#: 20910100A1
Last Changed By: J. HAWKINS #42401

154

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:32:37
Cal Check	0.101	09:32:56
Air Blank	0.000	09:33:25

Clare Conley
Operator's Signature