



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
5.31.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/11 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 4/28/11 and 5/26/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -


Sarah Cione

Sarah L. Cione

5/30/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30 day of MAY, 20 11.

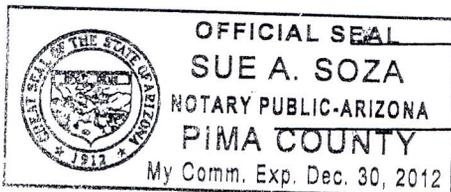

Sue A. Soza

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
5.31.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/26/11 TIME 856

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7sc

BH

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
05/26/2011 08:58:33

Standard Lot#: 20910100A1  
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:59:20
Diagnostic Test	Pass	08:59:48
Air Blank	0.000	09:00:15
0.100 Cal Check	0.101	09:00:34
Air Blank	0.000	09:01:03
Subject Test	0.000	09:01:42
Air Blank	0.000	09:02:09
Five Minute Wait	ABT*	
Air Blank	0.000	09:02:44

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7 sc

BJ

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
05/26/2011 09:02:49

Standard Lot#: 20910100A1  
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:03:13
Air Blank	IPS**	09:03:20

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7 sc

BY

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
05/26/2011 09:04:11

Standard Lot#: 20910100A1  
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:05:15
Diagnostic Test	Pass	09:05:42
Air Blank	0.000	09:06:09
0.100 Cal Check	0.100	09:06:28
Air Blank	0.000	09:06:58
Subject Test	0.033*	09:10:06
Air Blank	0.000	09:10:38
Air Blank	0.000	09:11:05
Subject Test	XXX**	09:11:24
Air Blank	0.000	09:11:58

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002264  
05/26/2011 09:16:49

5 of 7 SC

BSH

-----  
Standard Lot#: 20910100A1  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:17:16
Cal Check	0.101	09:17:35
Air Blank	0.000	09:18:04

*Sarah Cione*

Operator's Signature

6 of 7 sc

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TBY  
5.31.2011

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/26/11 TIME 856  
Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 75c

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002264  
05/26/2011 08:56:56

TS4

-----  
Standard Lot#: 20910100A1  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:57:23
Cal Check	0.101	08:57:42
Air Blank	0.000	08:58:11

Sarah Cione  
Operator's Signature