



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

FD4
7-26-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/21/11 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 6/23/11 and 7/21/11
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on
The instrument was returned to service on

Sarah Cione

Sarah L. Cione

7/25/11


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 25 day of July, 20 11.

A.M. Hernández

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

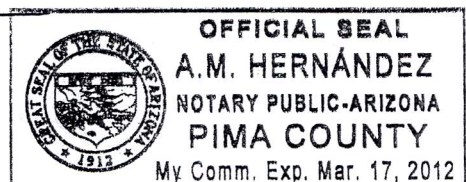

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
7-26-2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/21/11 TIME 905

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Sarah Cione

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

07/21/2011

09:06:35

Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:07:19
Diagnostic Test	Pass	09:07:47
Air Blank	0.000	09:08:14
0.100 Cal Check	0.098	09:08:33
Air Blank	0.000	09:09:02
Subject Test	0.000	09:09:50
Air Blank	0.000	09:10:17
Five Minute Wait	ABT*	
Air Blank	0.000	09:11:08

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

07/21/2011

09:15:19

Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:15:51
Air Blank	IPS**	09:15:58

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

07/21/2011

09:17:04

Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:18:02
Diagnostic Test	Pass	09:18:30
Air Blank	0.000	09:18:57
0.100 Cal Check	0.098	09:19:16
Air Blank	0.000	09:19:45
Subject Test	0.027*	09:22:54
Air Blank	0.000	09:23:25
Air Blank	0.000	09:23:52
Subject Test	XXX**	09:24:22
Air Blank	0.000	09:24:55

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

07/21/2011

09:25:48

Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:26:15
Cal Check	0.098	09:26:34
Air Blank	0.000	09:27:03

Sarah Cione
Operator's Signature

5 of 7_{sc}

PM

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4
7-26-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 7/21/11 TIME 905
Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

07/21/2011

09:05:00

Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:05:28
Cal Check	0.098	09:05:47
Air Blank	0.000	09:06:16



Operator's Signature

7 of 7 sc

Ph



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2
sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

BS
7-26-2011

Date: 6/24/11
QAS: HAWKINS 42401
Intoxilyzer Serial Number: 80-002264
New Lot Number: 04911100A1
New Expiration Date: 3/1/2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

[Signature] 42401

2 of 2sc


BH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
06/24/2011 23:26:42

Standard Lot#: 04911100A1
Last Changed By: J. HAWKINS #42401

QAS: J. HAWKINS #42401
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	23:27:09
Cal Check	0.099	23:27:28
Air Blank	0.000	23:27:58

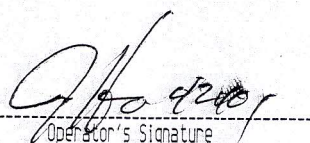
Operator's Signature