



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
11-7-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 9/15/11 and 10/14/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 20 11.

A.M. Hernandez

Notary Public

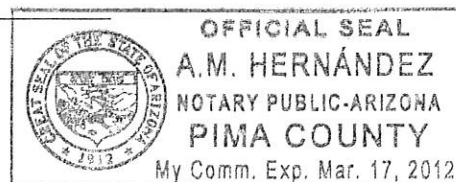

\*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7  
sc

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
11.7.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/11 TIME 1006

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7sc

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

10/14/2011

10:09:11

TS4

Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test             | g/210L | Time     |
|------------------|--------|----------|
| Air Blank        | 0.000  | 10:10:01 |
| Diagnostic Test  | Pass   | 10:10:28 |
| Air Blank        | 0.000  | 10:10:56 |
| 0.100 Cal Check  | 0.100  | 10:11:15 |
| Air Blank        | 0.000  | 10:11:44 |
| Subject Test     | 0.000  | 10:12:11 |
| Air Blank        | 0.000  | 10:12:38 |
| Five Minute Wait | ABT*   |          |
| Air Blank        | 0.000  | 10:13:19 |

\*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 75

TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
10/14/2011

10:14:20

Standard Lot#: 04911100A1  
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | RFI*   | 10:14:50 |
| Air Blank | IPS**  | 10:14:57 |

\*RFI Detect

\*\*Improper Sample

Not a Successfully Completed Test Sequence

4 of 7sc

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002264

10/14/2011

10:16:55

TS4

Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS

#42401

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:17:39 |
| Diagnostic Test | Pass   | 10:18:07 |
| Air Blank       | 0.000  | 10:18:34 |
| 0.100 Cal Check | 0.100  | 10:18:53 |
| Air Blank       | 0.000  | 10:19:22 |
| Subject Test    | 0.069* | 10:22:31 |
| Air Blank       | 0.000  | 10:23:02 |
| Air Blank       | 0.000  | 10:23:30 |
| Subject Test    | XXX**  | 10:24:14 |
| Air Blank       | 0.000  | 10:24:49 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7sc

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

10/14/2011

10:25:30

TSL

-----  
Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 10:25:58 |
| Cal Check | 0.098  | 10:26:16 |
| Air Blank | 0.000  | 10:26:46 |

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature

6 of 7sc

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TJL  
11.7.2011

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 10/14/11 TIME 1006  
Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7<sub>sc</sub>

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002264

10/14/2011

10:06:53

TS4

-----  
Standard Lot#: 04911100A1

Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 10:07:21 |
| Cal Check | 0.099  | 10:07:39 |
| Air Blank | 0.000  | 10:08:09 |



Operator's Signature