



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TSG  
 12.15.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-8-11 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 11-10-11 and 12-8-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.


Clare Conley  
 Clare L. Conley

12-8-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8<sup>th</sup>, day of December, 2011.

Ann L. Beecroft  
 Notary Public

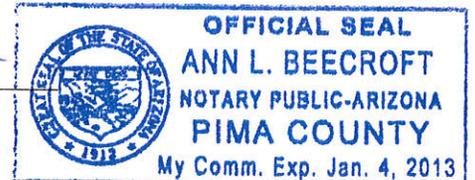

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

1068  
en  
TS4 12.15.2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 12-8-11 TIME 0717  
INTOXILYZER SERIAL # 82-002264 LOCATION ODW Day HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
12/08/2011

07:34:05

2568  
cm

---

Standard Lot#: 04911100A1  
Last Changed By: J. HAWKINS #42401

---

TS4

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	07:34:53
Diagnostic Test	Pass	07:35:20
Air Blank	0.000	07:35:48
0.100 Cal Check	0.100	07:36:06
Air Blank	0.000	07:36:36
Subject Test	0.000	07:37:12
Air Blank	ABT*	07:37:18
Air Blank	RFI**	07:37:25

\*Sequence Aborted  
\*\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
12/08/2011

07:37:55

308  
CW

---

Standard Lot#: 04911100A1  
Last Changed By: J. HAWKINS #42401

---

B4

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	IPS*	07:38:22
Air Blank	0.000	07:38:50

\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
12/08/2011 07:38:54

4578  
em  
B4

Standard Lot#: 04911100A1  
Last Changed By: J. HAWKINS #42401

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST,SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:39:37
Diagnostic Test	Pass	07:40:05
Air Blank	0.000	07:40:32
0.100 Cal Check	0.100	07:40:51
Air Blank	0.000	07:41:20
Subject Test	0.221*	07:44:29
Air Blank	0.000	07:45:03
Air Blank	0.000	07:45:30
Subject Test	RGE**	07:46:02
Air Blank	0.000	07:46:37

\*Deficient Sample  
\*\*Range Exceeded

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002264  
12/08/2011 07:46:42

508  
OK

Standard Lot#: 04911100A1  
Last Changed By: J. HAWKINS #42401

BS4

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:47:25
Diagnostic Test	Pass	07:47:52
Air Blank	0.000	07:48:19
0.100 Cal Check	0.100	07:48:38
Air Blank	0.000	07:49:08
Subject Test	XXX*	07:49:22
Air Blank	0.000	07:49:56

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002264  
12/08/2011 07:50:14

6.68  
Qu  
T39

-----  
Standard Lot#: 04911100A1  
Last Changed By: J. HAWKINS #42401  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:50:41
Cal Check	0.099	07:51:00
Air Blank	0.000	07:51:30

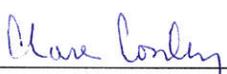
  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

cn  
~~8~~  
7568 cn

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
12-15-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 12-8-11 TIME 0717  
Intoxilyzer Serial # 80-802264 Location ODW Jail HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002264  
12/08/2011 07:17:21

8068  
ew

-----  
Standard Lot#: 04911100A1  
Last Changed By: J. HAWKINS #42401  
-----

TS4

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:17:49
Cal Check	0.101	07:18:08
Air Blank	0.000	07:18:37

*Clare Conley*

\_\_\_\_\_  
Operator's Signature