



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
3-20-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/1/12 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 2/2/12 and 3/1/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 3/1/12. The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


3/9/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th day of March, 2012.

Elma Escarcega

Notary Public

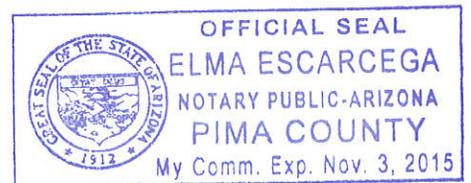

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 10_{sc}

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSG
3.20.2012

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/1/12 TIME 1053

INTOXILYZER SERIAL # 80-002264 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

2 of 10_{sc}
TBY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
03/01/2012 10:56:06

Standard Lot#: 04911100A1
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 10:56:58 |
| Diagnostic Test | Pass | 10:57:26 |
| Air Blank | 0.000 | 10:57:53 |
| 0.100 Cal Check | 0.101 | 10:58:12 |
| Air Blank | 0.000 | 10:58:41 |
| Subject Test | 0.000 | 10:59:40 |
| Air Blank | 0.000 | 11:00:07 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 11:01:23 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 10_{sc}

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
03/01/2012 11:06:14

Standard Lot#: 04911100A1
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:07:04 |
| Diagnostic Test | Pass | 11:07:32 |
| Air Blank | 0.000 | 11:07:59 |
| 0.100 Cal Check | RFI* | 11:08:10 |
| Air Blank | IPS** | 11:08:19 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 10_{sc}
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002264
03/01/2012 11:08:25

Standard Lot#: 04911100A1
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:09:08 |
| Diagnostic Test | Pass | 11:09:35 |
| Air Blank | 0.000 | 11:10:03 |
| 0.100 Cal Check | 0.101 | 11:10:21 |
| Air Blank | 0.000 | 11:10:51 |
| Subject Test | 0.063* | 11:13:59 |
| Air Blank | 0.000 | 11:14:31 |
| Air Blank | 0.000 | 11:14:58 |
| Subject Test | XXX** | 11:15:28 |
| Air Blank | 0.000 | 11:16:04 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 10_{sc}


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/01/2012 11:16:24

TS4

Standard Lot#: 04911100A1
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 11:16:52 |
| Cal Check | 0.100 | 11:17:11 |
| Air Blank | 0.000 | 11:17:40 |

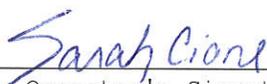

Operator's Signature

EXHIBIT G-2

6 of 10
sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

TS4
3-20-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/1/12 TIME 1053
Intoxilyzer Serial # 80-002264 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓)✓ 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/01/2012 10:53:05

7 of 10 sc
BH

Standard Lot#: 04911100A1
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 10:53:32 |
| Cal Check | INT* | 10:53:51 |
| Air Blank | 0.000 | 10:54:22 |

*Interferent Detect



Operator's Signature

8 of 10_{sc}


TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002264
03/01/2012 10:54:30

Standard Lot#: 04911100A1
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 10:54:58 |
| Cal Check | 0.101 | 10:55:16 |
| Air Blank | 0.000 | 10:55:46 |

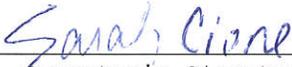
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

9 of 10_{sc}

INTOXILYZER SERVICE RECORD

T34
3.20.2012

Date: 3/1/12

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002264

Reason For Service: preventative maintenance

Date Out of Service: 3/1/12

Calibration Check completed prior to removal from service

Date: 3/1/12

SQAP completed prior to removal from service

Date: 3/1/12

Corrective Action: Sent to CMI, Inc.



Service Evaluation Form

10 of 10 sc

TJH

This form **MUST** be completed and enclosed with instrument to be serviced.
Failure to complete and return this form may cause delays in service.

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)
 Name Terry Gallegos Phone: (520) 837-7677
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:
Finance / Accounts Payable Tucson Police Department Crime Lab
Po Box 27450 1306 W. Miracle Mile
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-002264 Instrument Model: I8000

4. Detailed Description of Problem:
preventative maintenance + software update

*****Hazardous Material Warning! - DO NOT return gas cylinders with instrument!*****

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:
Terry Gallegos Toxicology Coordinator
 Name (Please Print) Title
Terry Gallegos 3.20.2012
 Signature Date

Ship to:
 CMI, Inc.
 Attn: Service Dept.
 316 East Ninth Street
 Owensboro, KY 42303

No, please send estimate before repairs are made.
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.