



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
6/18/12

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6.11.2012 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of — and — the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on 6.11.2012.

Terry S. Gallegos

Terry S. Gallegos


6.13.2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13th day of June, 2012.

Elma Escarcega

Notary Public

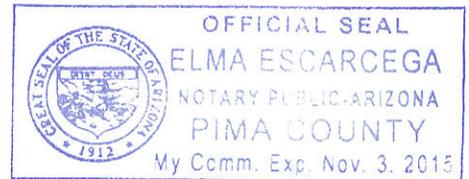

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

10 of 14
TS4
SC
6/18/12

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-002264

Date Returned to Service: 6.11.2012

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

20f14
TS4
SC
6/18/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6.11.2012 TIME 1649
INTOXILYZER SERIAL # 80-002264 LOCATION CrimeLab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS _____

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
06/11/2012 16:51:52

3 of 14
TJ
SC

Standard Lot#: 04911100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:52:48
Diagnostic Test	Pass	16:53:23
Air Blank	0.000	16:53:51
0.100 Cal Check	0.098	16:54:12
Air Blank	0.000	16:54:41
Subject Test	0.000	16:55:00
Air Blank	0.000	16:55:28
Wait	ABT*	
Air Blank	0.000	16:56:02

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
06/11/2012 16:56:55

40F14
TS4

Standard Lot#: 04911100A1
Last Changed By: T. GALLEGOS #28847

5C

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:57:46
Diagnostic Test	Pass	16:58:21
Air Blank	0.000	16:58:49
0.100 Cal Check	0.098	16:59:09
Air Blank	RFI*	16:59:16
Air Blank	IPS**	16:59:23

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
06/11/2012

16:59:55

50f14
B4

Standard Lot#: 04911100A1
Last Changed By: T. GALLEGOS #28847

5C

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:00:46
Diagnostic Test	Pass	17:01:21
Air Blank	0.000	17:01:48
0.100 Cal Check	0.099	17:02:09
Air Blank	0.000	17:02:39
Subject Test	DEF*	17:05:46
Air Blank	0.000	17:06:14
Air Blank	0.000	17:06:41
Subject Test	XXX**	17:07:07
Air Blank	0.000	17:07:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
06/11/2012 17:08:01

60 of 14
TS4
SC

Standard Lot#: 04911100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:08:32
0.100 Cal Check	0.097	17:08:52
Air Blank	0.000	17:09:22

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

70514
TS4
SC

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6.11.2012 TIME 1649
Intoxilyzer Serial # 80-002264 Location CrimeLab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
06/11/2012 16:49:58

Standard Lot#: 04911100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:50:30
0.100 Cal Check	0.099	16:50:50
Air Blank	0.000	16:51:20

Preventative Maintenance
Performed -- Timer Reset

80514
TS4
SC



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 6.11.2012 9 of 14
TS4 SC
6/18/12
 QAS: Terry Gallegos
 Intoxilyzer Serial Number: 80-002264
 New Lot Number: 04911100A1
 New Expiration Date: 3.1.2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press
 ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove ~~and sign printed record~~ ^{TS4}. Attach the record to the completed checklist.
 used same test for G-2.
 TS4

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
06/11/2012 16:49:58

Standard Lot#: 04911100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test g/210L Time

Air Blank 0.000 16:50:30
0.100 Cal Check 0.099 16:50:50
Air Blank 0.000 16:51:20

Preventative Maintenance
Performed -- Timer Reset

100814
TSh
SC

12 of 14
TS4

SC

Certificate of Calibration

Intoxilyzer® 8000

This is to certify that the calibration of **Intoxilyzer**® serial number 80-002264, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.

Date May 3, 2012 Signed D. DeRohrer

Technician

EMI INC.

316 East 9th Street
Owensboro, KY 42303
USA

Part No. 650519 CMI 7/9/09

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



Service Evaluation Form

130814
DJ

This form **MUST** be completed and enclosed with instrument to be serviced.
Failure to complete and return this form may cause delays in service.

SC

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)
Name Terry Gallegos Phone: (520) 837-7677
Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:
Finance / Accounts Payable Tucson Police Department Crime Lab
Po Box 27450 1306 W. Miracle Mile
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-002264 Instrument Model: Intox 8000

4. Detailed Description of Problem: Instrument just returned from service. Not registering any pressure in gas tank. Replaced regulator, still registers 0 psi. → see attached

Hazardous Material Warning! - DO NOT return gas cylinders with instrument!

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:
Terry Gallegos Toxicology Coordinator
Name (Please Print) Title
Terry Gallegos 5.14.2012
Signature Date

Ship to:
CMI, Inc.
Attn: Service Dept.
316 East Ninth Street
Owensboro, KY 42303

No, please send estimate before repairs are made.
Note: An estimate will be faxed before performing any repairs and may cause delays in service.
An evaluation fee (\$79.00 or actual costs) will apply to estimates

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

