



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

By  
7-24-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/20/12 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 6/20/12 and 7/20/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 7/20/12.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


7/23/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 23<sup>rd</sup>, day of July, 2012.

Elma Escarcega

Notary Public

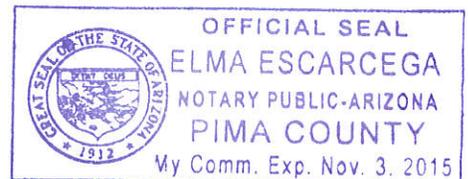

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 9 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

BA  
7-24-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/20/12 TIME 1616

INTOXILYZER SERIAL # 80-002264 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0.000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 95c  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
07/20/2012 16:18:01

---

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:18:43
Diagnostic Test	Pass	16:19:18
Air Blank	0.000	16:19:46
0.100 Cal Check	0.098	16:20:06
Air Blank	0.000	16:20:36
Subject Test	0.000	16:21:13
Air Blank	0.000	16:21:40
Wait	ABT*	
Air Blank	0.000	16:22:14

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 9  
sc  
ASH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
07/20/2012 16:22:24

---

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:23:06
Diagnostic Test	Pass	16:23:41
Air Blank	0.000	16:24:09
0.100 Cal Check	0.097	16:24:29
Air Blank	0.000	16:24:59
Subject Test	RFI*	16:25:29
Air Blank	IPS**	16:25:37

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
07/20/2012

16:25:42

4 of 9 sc  
BSH

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:26:24
Diagnostic Test	Pass	16:26:59
Air Blank	0.000	16:27:26
0.100 Cal Check	0.097	16:27:47
Air Blank	0.000	16:28:16
Subject Test	0.079*	16:31:24
Air Blank	0.000	16:31:56
Air Blank	0.000	16:32:24
Subject Test	XXX**	16:32:40
Air Blank	0.000	16:33:15

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
07/20/2012 16:33:26

5 of 9 SC  
TS4

-----  
Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:33:56
0.100 Cal Check	0.097	16:34:16
Air Blank	0.000	16:34:46

Preventative Maintenance NOT  
Performed -- Timer NOT Reset

EXHIBIT G-2

6 of 9sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
7-24-2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 7/20/12 TIME 1616  
Intoxilyzer Serial # 80-002264 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
07/20/2012 16:16:07

7 of 9 SC

154

-----  
Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:16:37
0.100 Cal Check	0.097	16:16:57
Air Blank	0.000	16:17:27

Preventative Maintenance  
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD

TS4  
7.24.2012

Date: 7/20/12

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002264

Reason For Service: voltage/current failures when plugged into DC, periodic "Contact QAS- tank pressure below minimum"

Date Out of Service: 7/20/12

Calibration Check completed prior to removal from service

Date: 7/20/12

SQAP completed prior to removal from service

Date: 7/20/12

Corrective Action:  Sent to CMI, Inc.

data downloaded



# Service Evaluation Form

9 of 9 SC  
TS4

This form MUST be completed and enclosed with instrument to be serviced.  
Failure to complete and return this form may cause delays in service.

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)  
Name Terry Gallegos Phone: (520) 837-7677  
Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Finance / Accounts Payable  
Po Box 27450  
Tucson, AZ 85726-7450  
Ship to Address: Tucson Police Department Crime Lab  
1306 W. Miracle Mile  
Tucson, AZ 85705

3. Serial Number: 80-002264 Instrument Model: I8000

4. Detailed Description of Problem:  
voltage/current failures when plugged into DC, periodic "Contact QAS-  
tank pressure below minimum"

\*\*\*Hazardous Material Warning! - DO NOT return gas cylinders with instrument!\*\*\*

I Authorize Repairs Up To:  All  \$250  \$500  \$750  Other \$ \_\_\_\_\_  
Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator  
Name (Please Print) Title  
Terry Gallegos 7-24-2012  
Signature Date

Ship to:  
CMI, Inc.  
Attn: Service Dept.  
316 East Ninth Street  
Owensboro, KY 42303

No, please send estimate before repairs are made.  
Note: An estimate will be faxed before performing any repairs and may cause delays in service.  
An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.