



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
9.20.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/14/12 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of — and — the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on 9/14/12.

Sarah Cione

Sarah L. Cione


9/17/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of September, 2012.

Elma Escarcega

Notary Public

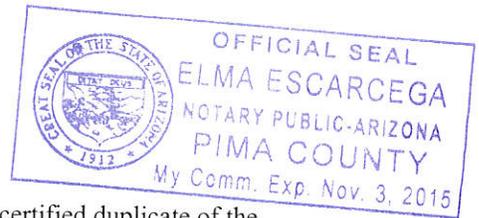

*Includes documents specified in A.A.C. R13-10-104.A.5



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

1 of 11 sc

TS4
9-20-2012

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002264

Date Returned to Service: 9/14/12

2 of 11 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
9.20.2012

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/14/12 TIME 1407

INTOXILYZER SERIAL # 80-002264 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 11 sc

134

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
09/14/2012 14:09:27

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:10:10
Diagnostic Test	Pass	14:10:45
Air Blank	0.000	14:11:12
0.100 Cal Check	0.100	14:11:33
Air Blank	0.000	14:12:03
Subject Test	0.000	14:12:17
Air Blank	0.000	14:12:44
Wait	ABT*	
Air Blank	0.000	14:13:22

*Sequence Aborted

Not a Successfully Completed Test Sequence

4 of 11 sc
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
09/14/2012 14:13:28

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:14:11
Diagnostic Test	Pass	14:14:46
Air Blank	0.000	14:15:13
0.100 Cal Check	0.100	14:15:34
Air Blank	0.000	14:16:04
Subject Test	RFI*	14:16:20
Air Blank	IPS**	14:16:28

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

5 of 11 sc

TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
09/14/2012 14:16:34

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:17:17
Diagnostic Test	Pass	14:17:52
Air Blank	0.000	14:18:19
0.100 Cal Check	0.100	14:18:40
Air Blank	0.000	14:19:10
Subject Test	0.045*	14:22:17
Air Blank	0.000	14:22:50
Air Blank	0.000	14:23:17
Subject Test	XXX**	14:23:38
Air Blank	0.000	14:24:12

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

6 of 11 sc

B4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
09/14/2012 14:24:24

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:24:52
0.100 Cal Check	0.099	14:25:13
Air Blank	0.000	14:25:43

Preventative Maintenance
Performed -- Timer Reset

7 of 11 SC
TSH
9.20.2012

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/14/12 TIME 1407
Intoxilyzer Serial # 80-002264 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
09/14/2012 14:07:49

8 of 11sc
B4

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:08:18
0.100 Cal Check	0.100	14:08:39
Air Blank	0.000	14:09:09

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

9 of 11 sc
TS4

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 9/14/12
 QAS: Sarah Cione
 Intoxilyzer Serial Number: 80-002264
 New Lot Number: 13712100A2
 New Expiration Date: 7/1/2014

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove ~~and sign~~ printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

10 of 11 SC
TSG

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
09/14/2012 14:05:17

Standard Lot#: 13712100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:05:46
0.100 Cal Check	0.100	14:06:07
Air Blank	0.000	14:06:36

Preventative Maintenance
Performed -- Timer Reset

