



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

SC  
 9/16/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-22-2013 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 7-25-2013 and 8-22-2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on —.

Terry Gallegos

Terry S. Gallegos


9-10-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 10, day of September, 2013.

A. M. Hernandez

Notary Public

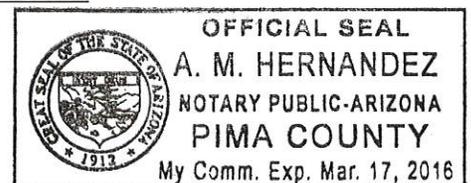

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TSg  
sc  
9/16/13

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 8.22.2013 TIME 1042  
INTOXILYZER SERIAL # 80-002264 LOCATION ODM

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
08/22/2013

10:44:36

2077  
Dh

Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616

Sc

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:45:26
Diagnostic Test	Pass	10:46:01
Air Blank	0.000	10:46:29
0.100 Cal Check	0.100	10:46:50
Air Blank	0.000	10:47:20
Subject Test	0.000	10:47:53
Air Blank	0.000	10:48:20
Wait	ABT*	
Air Blank	0.000	10:48:55

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
08/22/2013

10:49:14

30ft  
D4

Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616

SC

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:50:16
Diagnostic Test	Pass	10:50:51
Air Blank	0.000	10:51:18
0.100 Cal Check	RFI*	10:51:32
Air Blank	IPS**	10:51:42

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
08/22/2013

10:51:48

4 of 7  
TSg

Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:52:54
Diagnostic Test	Pass	10:53:29
Air Blank	0.000	10:53:57
0.100 Cal Check	0.100	10:54:17
Air Blank	0.000	10:54:48
Subject Test	0.155*	10:57:55
Air Blank	0.000	10:58:29
Air Blank	0.000	10:58:56
Subject Test	XXX**	10:59:23
Air Blank	0.000	10:59:59

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
08/22/2013 11:00:16

Soft  
TS4

-----  
Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616  
-----

SC

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:00:47
0.100 Cal Check	0.098	11:01:08
Air Blank	0.000	11:01:38

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TS4  
SC  
9/16/13

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 8-22-2013 TIME 1042  
Intoxilyzer Serial # 80-002264 Location ODM

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
08/22/2013 10:42:29

7 of 7  
TS4

SC

-----  
Standard Lot#: 13712100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	10:42:56
0.100 Cal Check	0.100	10:43:17
Air Blank	0.000	10:43:47