



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

SC
5/19/14

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.1.2014 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 4.3.2014 and 5.1.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Terry Gallegos

Terry S. Gallegos


5.15.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 15th, day of May, 2014.

Edma Escarcega

Notary Public

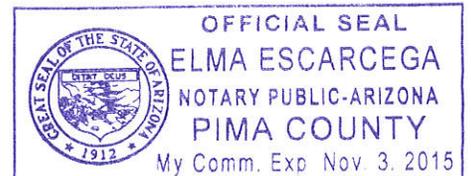

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10f9
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
5/19/14

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 5.1.2014 TIME 1042

INTOXILYZER SERIAL # 80-002264 LOCATION ODM

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
05/01/2014

10:44:15

20f9
TS4
SC

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:45:09
Diagnostic Test	Pass	10:45:44
Air Blank	0.000	10:46:12
0.100 Cal Check	0.101	10:46:33
Air Blank	0.000	10:47:03
Subject Test	RFI*	10:47:23
Air Blank	RFI*	10:47:31

*RFI Detect

Due to cell phone TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
05/01/2014

10:47:42

30f9
TS4
SC

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:48:18
Air Blank	0.000	10:48:47

*RFI Detect

Turned Cell phone off after this TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
05/01/2014

10:49:15

40f9
TSg
sc

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:50:12
Diagnostic Test	Pass	10:50:47
Air Blank	0.000	10:51:14
0.100 Cal Check	0.101	10:51:35
Air Blank	0.000	10:52:05
Subject Test	0.000	10:52:30
Air Blank	0.000	10:52:57
Wait	ABT*	
Air Blank	0.000	10:53:32

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
05/01/2014 10:53:38

50f9
TS4
SC

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:54:27
Diagnostic Test	Pass	10:55:02
Air Blank	0.000	10:55:29
0.100 Cal Check	RFI*	10:55:43
Air Blank	IPS**	10:55:51

*RFI Detect
**Improper Sample

TPD Radio TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
05/01/2014

10:55:56

60f9
TSg

SC

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:56:45
Diagnostic Test	Pass	10:57:20
Air Blank	0.000	10:57:48
0.100 Cal Check	0.101	10:58:08
Air Blank	0.000	10:58:38
Subject Test	0.271*	11:01:46
Air Blank	0.000	11:02:20
Air Blank	0.000	11:02:48
Subject Test	XXX**	11:03:08
Air Blank	0.000	11:03:45

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
05/01/2014 11:04:01

70F9
TS4

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

SC

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:04:30
0.100 Cal Check	0.099	11:04:51
Air Blank	0.000	11:05:21

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

80f9
TS4
SC
5/19/14

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5-1-2014 TIME 1042
Intoxilyzer Serial # 80-002264 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
05/01/2014 10:42:34

90f9
TS4
SC

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	10:43:01
0.100 Cal Check	0.101	10:43:22
Air Blank	0.000	10:43:52