



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
1-21-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1.7.2015 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 12.11.2014 and 1.7.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


1/8/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 8th, day of January, 2015.

Elma Escarcega

Notary Public

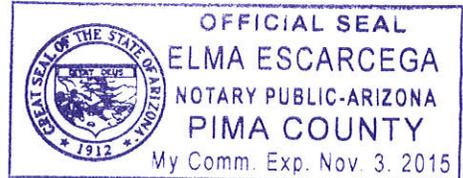

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
WLL

By
1.21.2015

QA SPECIALIST WILL LOOMIS 53545 AGENCY Tucson Police Department
DATE 1/7/15 TIME 08:33
INTOXILYZER SERIAL # 80-002264 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS RFI DETECTED DURING CALIBRATION CHECK;
CALIBRATION CHECK REPEATED

SIGNATURE Will Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
01/07/2015

08:37:02

2 of 8
WLL
TS4

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1980
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:38:13
Diagnostic Test	Pass	08:38:48
Air Blank	0.000	08:39:16
0.100 Cal Check	0.100	08:39:36
Air Blank	0.000	08:40:06
Subject Test	0.000	08:40:27
Air Blank	0.000	08:40:54
Wait	ABT*	
Air Blank	0.000	08:41:29

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
01/07/2015

08:41:52

3 of 8
WLL
TSH

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:42:59
Diagnostic Test	Pass	08:43:34
Air Blank	0.000	08:44:02
0.100 Cal Check	RFI*	08:44:15
Air Blank	IPS**	08:44:24

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
01/07/2015

08:47:01

4 of 8
WLL
TSY

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:47:59
Diagnostic Test	Pass	08:48:34
Air Blank	0.000	08:49:01
0.100 Cal Check	0.100	08:49:22
Air Blank	0.000	08:49:52
Subject Test	0.162*	08:52:59
Air Blank	0.000	08:53:33
Air Blank	0.000	08:54:00
Subject Test	XXX**	08:54:16
Air Blank	0.000	08:54:50

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
01/07/2015

08:44:51

5 of 8 WLL
BM

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:45:53
Diagnostic Test	RFI*	08:46:17
Air Blank	0.000	08:46:46

* RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
01/07/2015 08:55:58

6 of 8
WLL

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

TJH

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:56:30
0.100 Cal Check	0.099	08:56:50
Air Blank	0.000	08:57:20

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8
WLU
134
1.21.2015

QA Specialist WILL LOOMIS 53545 AGENCY Tucson Police Department
DATE 1/7/15 TIME 08:33
Intoxilyzer Serial # 80-002264 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Core Version: 8105.48
Serial Number: 80-002264
01/07/2015 08:33:58

8 of 8
WLL
TS4

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:34:31
0.100 Cal Check	0.100	08:34:52
Air Blank	0.000	08:35:22
End of test		

**Preventive Maintenance NOT
Performed -- Timer NOT Reset**

Reprinted From Database

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.