



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
 4.10.2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/2/2015 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 3/5/2015 and 4/2/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.

William L Loomis

William L. Loomis


4/2/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 2nd, day of April, 2015.

Elma Escarcega

Notary Public

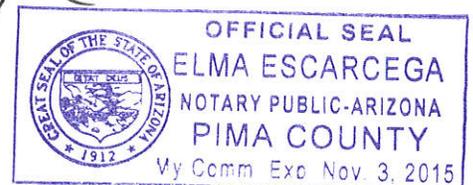

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
WLL
TS 4.10.2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 4/2/2015 TIME 08:01

INTOXILYZER SERIAL # 80-002264 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS RFI DETECTED JUST AFTER DEFICIENT SAMPLE.
TEST REPEATED FOR DEFICIENT SAMPLE.

SIGNATURE Will Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
04/02/2015

08:03:35

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST

DOB: 11/11/1987

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:04:28
Diagnostic Test	Pass	08:05:03
Air Blank	0.000	08:05:31
0.100 Cal Check	0.100	08:05:51
Air Blank	0.000	08:06:21
Subject Test	0.000	08:06:39
Air Blank	0.000	08:07:06
Wait	ABT*	
Air Blank	0.000	08:07:43

*Sequence Aborted

2 of 8
WLL
BY

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
04/02/2015

08:07:49

3 of 8
WU
154

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:08:38
Diagnostic Test	Pass	08:09:13
Air Blank	0.000	08:09:41
0.100 Cal Check	RFI*	08:09:55
Air Blank	IPS**	08:10:03

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
04/02/2015

08:10:08

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:11:02
Diagnostic Test	Pass	08:11:37
Air Blank	0.000	08:12:05
0.100 Cal Check	0.100	08:12:25
Air Blank	0.000	08:12:55
Subject Test	0.400*	08:16:05
Air Blank	RFI**	08:16:11
Air Blank	0.000	08:16:41

*Deficient Sample
**RFI Detect

4 of 8
WLL
B4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
04/02/2015

08:16:47

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:17:39
Diagnostic Test	Pass	08:18:14
Air Blank	0.000	08:18:41
0.100 Cal Check	0.099	08:19:02
Air Blank	0.000	08:19:32
Subject Test	0.222*	08:22:40
Air Blank	0.000	08:23:14
Air Blank	0.000	08:23:42
Subject Test	XXX**	08:24:05
Air Blank	0.000	08:24:40

*Deficient Sample

**Invalid Sample

Begin new deprivation period

5 of 8
WU
TS4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
04/02/2015 08:25:08

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:25:41
0.100 Cal Check	0.099	08:26:02
Air Blank	0.000	08:26:32

Preventative Maintenance
Performed -- Timer Reset

6 of 8
WU
Bh

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 4/2/2015 TIME 08:01
Intoxilyzer Serial # 80-002264 Location ODM

7 of 8
WLL
TS4
4.10.2015

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002264
Core Version: 8105.48
04/02/2015 08:01:58

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:02:31
0.100 Cal Check	0.101	08:02:51
Air Blank	0.000	08:03:21

Preventative Maintenance
Performed -- Timer Reset

8 of 8
WLL
TSL