



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/28/2015 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 4/30/2015 and 5/28/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on
The instrument was returned to service on

Handwritten signature of William L. Loomis

William L. Loomis


6/8/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 8th day of June, 2015.

Handwritten signature of Notary Public

Notary Public

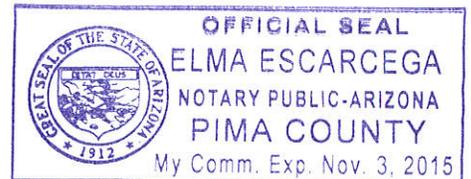

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8  
WLL

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 5.28.15 TIME 08:10

INTOXILYZER SERIAL # 80-002264 LOCATION ODM

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(  ) 1. Clock time check.

(  ) 2. Date check.

OPERATIONAL TESTS

(  ) 1. Alcohol-free subject test result 0. 000 AC.

(  ) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(  ) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(  ) 4. Standard Calibration Check standard 0. 100 AC

Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS RFI DETECTED DURING STANDARD CALIBRATION  
CHECK, CALIBRATION CHECK REPEATED

SIGNATURE Will Loomis

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
05/28/2015

08:13:14

2 of 8  
WLL

---

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

---

OAS: W. LOOMIS #53545  
TUCSON PD

---

Operator: W. LOOMIS #53545  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:14:06
Diagnostic Test	Pass	08:14:41
Air Blank	0.000	08:15:09
0.100 Cal Check	0.101	08:15:30
Air Blank	0.000	08:16:00
Subject Test	0.000	08:16:18
Air Blank	0.000	08:16:46
Wait	ABT*	
Air Blank	0.000	08:17:20

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
05/28/2015

08:18:16

3 of 8  
WLL

---

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: W. LOOMIS #53545  
TUCSON PD

---

Operator: W. LOOMIS #53545  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:19:07
Diagnostic Test	Pass	08:19:42
Air Blank	0.000	08:20:10
0.100 Cal Check	RFI*	08:20:22
Air Blank	IPS**	08:20:31

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
05/28/2015

08:20:38

4 of 8  
WLL

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:21:30
Diagnostic Test	Pass	08:22:05
Air Blank	0.000	08:22:33
0.100 Cal Check	0.100	08:22:53
Air Blank	0.000	08:23:23
Subject Test	0.075*	08:26:31
Air Blank	0.000	08:27:05
Air Blank	0.000	08:27:32
Subject Test	XXX**	08:27:49
Air Blank	0.000	08:28:25

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
05/28/2015 08:28:43

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:29:14
0.100 Cal Check	0.099	08:29:35
Air Blank	RFI*	08:29:49
Air Blank	0.000	08:30:17

\*RFI Detect

5 of 8  
WLL

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
05/28/2015 08:30:21

6 of 8  
WLL

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:30:51
0.100 Cal Check	0.099	08:31:12
Air Blank	0.000	08:31:42

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8  
will

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 5.28.2015 TIME 08:10  
Intoxilyzer Serial # 80-002264 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
05/28/2015

08:10:41

8 of 8  
WLU

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:11:13
0.100 Cal Check	0.100	08:11:33
Air Blank	0.000	08:12:03

Preventative Maintenance  
Performed -- Timer Reset