



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
7-8-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/25/2015 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 5/28/2015 and 6/25/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on [blank]
The instrument was returned to service on [blank]

William L. Loomis

William L. Loomis


7/6/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 6th, day of July, 2015.

Elma Escarcega

Notary Public

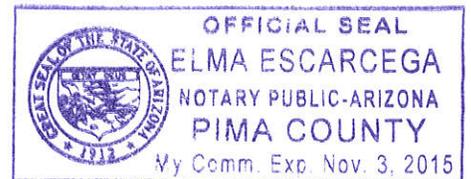

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

*1 of 7  
WLL  
TSH 7.8.2015*

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 6/25/2015 TIME 08:48  
INTOXILYZER SERIAL # 80-002264 LOCATION ODM

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
06/25/2015

08:52:08

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:53:00
Diagnostic Test	Pass	08:53:35
Air Blank	0.000	08:54:02
0.100 Cal Check	0.099	08:54:23
Air Blank	0.000	08:54:53
Subject Test	0.000	08:55:11
Air Blank	0.000	08:55:39
Wait	ABT*	
Air Blank	0.000	08:56:20

\*Sequence Aborted

2 of 7  
T84  
WLV

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
06/25/2015 08:56:33

3 of 7  
wlc  
TS4

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:57:24
Diagnostic Test	Pass	08:57:59
Air Blank	0.000	08:58:26
0.100 Cal Check	RFI*	08:58:40
Air Blank	IPS**	08:58:49

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
06/25/2015 08:58:54

4 of 7  
TS4  
WU

---

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: W. LOOMIS #53545  
TUCSON PD

---

Operator: W. LOOMIS #53545  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:59:46
Diagnostic Test	Pass	09:00:21
Air Blank	0.000	09:00:48
0.100 Cal Check	0.098	09:01:09
Air Blank	0.000	09:01:39
Subject Test	0.165*	09:04:46
Air Blank	0.000	09:05:20
Air Blank	0.000	09:05:48
Subject Test	XXX**	09:06:37
Air Blank	0.000	09:07:13

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
06/25/2015 09:07:48

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:08:18
0.100 Cal Check	0.098	09:08:39
Air Blank	0.000	09:09:09

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
WLL  
TJH

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLU  
By 7-8-2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 6/25/2015 TIME 08:48  
Intoxilyzer Serial # 80-002264 Location ODM

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
06/25/2015 08:48:36

7 of 7  
WLL  
TSh

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:49:06
0.100 Cal Check	0.099	08:49:27
Air Blank	0.000	08:49:57

Preventative Maintenance  
Performed -- Timer Reset