



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002264 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
9-22-2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-16-16 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 8-18-16 and 9-16-16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

[Handwritten signature of Amy L. Autret]

Amy L. Autret


9-19-16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 19th day of September, 2014.

[Handwritten signature of Elma Escarcega]

Notary Public

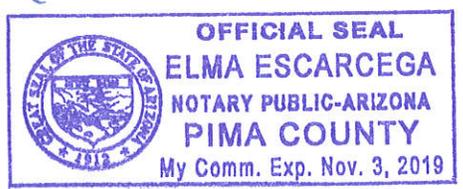

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TSK 9-22-2016

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 9-16-16 TIME 1240

INTOXILYZER SERIAL # 80-002264 LOCATION Laboratory

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS


- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

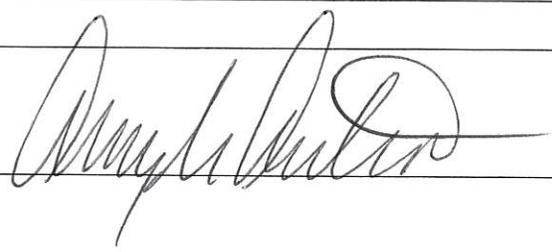

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
09/16/2016

12:42:27

2 of 7  
MS  
BY

---

Standard Lot#: 06616100A1  
Last Changed By: A. AUTRET #53546

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L           | Time     |
|-----------------|------------------|----------|
| Air Blank       | 0.000            | 12:43:18 |
| Diagnostic Test | Pass             | 12:43:53 |
| Air Blank       | 0.000            | 12:44:20 |
| 0.100 Cal Check | 0.097            | 12:44:41 |
| Air Blank       | 0.000            | 12:45:11 |
| Subject Test    | 0.000            | 12:45:27 |
| Air Blank       | 0.000            | 12:45:54 |
| Wait            | ABT <sup>^</sup> |          |
| Air Blank       | 0.000            | 12:46:30 |

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
09/16/2016

12:46:53

3 of 7  
ML  
TS4

---

Standard Lot#: 06616100A1  
Last Changed By: A. AUTRET #53546

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:47:44 |
| Diagnostic Test | Pass   | 12:48:19 |
| Air Blank       | 0.000  | 12:48:46 |
| 0.100 Cal Check | RFI*   | 12:49:00 |
| Air Blank       | IPS**  | 12:49:09 |

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
09/16/2016

12:49:15

4 of 7  
JLH  
TSH

---

Standard Lot#: 06616100A1  
Last Changed By: A. AUTRET #53546

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:50:07 |
| Diagnostic Test | Pass   | 12:50:42 |
| Air Blank       | 0.000  | 12:51:09 |
| 0.100 Cal Check | 0.098  | 12:51:30 |
| Air Blank       | 0.000  | 12:52:00 |
| Subject Test    | 0.228* | 12:55:07 |
| Air Blank       | 0.000  | 12:55:43 |
| Air Blank       | 0.000  | 12:56:11 |
| Subject Test    | XXX**  | 12:56:34 |
| Air Blank       | 0.000  | 12:57:10 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
09/16/2016 12:57:22

5067  
JL  
T84

-----  
Standard Lot#: 06616100A1  
Last Changed By: A. AUTRET #53546  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:57:50 |
| 0.100 Cal Check | 0.096  | 12:58:11 |
| Air Blank       | 0.000  | 12:58:42 |

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TS4  
9.22.2016

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department  
DATE 9-16-16 TIME 1240  
Intoxilyzer Serial # 80-002264 Location Laboratory

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002264  
Core Version: 8105.48  
09/16/2016 12:40:58

2 of 7  
[Handwritten signature]  
TJH

-----  
Standard Lot#: 06616100A1  
Last Changed By: A. AUTRET #53546  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:41:27 |
| 0.100 Cal Check | 0.097  | 12:41:47 |
| Air Blank       | 0.000  | 12:42:17 |

Preventative Maintenance  
Performed -- Timer Reset