



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

je

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-26-2006 Cal 1 0.099 AC Cal 2 0.100 AC


The instrument was placed into service on 12.26.2006.

Terry Gallegos
Terry S. Gallegos

1.4.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

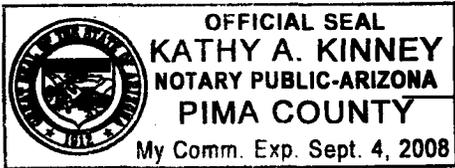

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

10 of 7
B4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

jk

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 12-26-2006 TIME 1544

INTOXILYZER SERIAL # 80-002265 LOCATION Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓
Deficient sample printed. ✓
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7
TB

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
12/26/2006 15:47:13

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:48:06
Diagnostic Test	Pass	15:48:33
Air Blank	0.000	15:49:01
0.100 Cal Check	0.100	15:49:19
Air Blank	0.000	15:49:49
Subject Test	0.000	15:50:09
Air Blank	0.000	15:50:37
Five Minute Wait	ABT*	
Air Blank	0.000	15:51:17

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

30F7
TSh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
12/26/2006

jk

15:51:27

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:52:19
Diagnostic Test	Pass	15:52:47
Air Blank	0.000	15:53:14
0.100 Cal Check	0.100	15:53:33
Air Blank	IPS*	15:53:41
Air Blank	RFI**	15:54:00

*Improper Sample
**RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT *ju*
Serial Number : 80-002265
12/26/2006 15:56:42

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:57:33
Diagnostic Test	Pass	15:58:00
Air Blank	0.000	15:58:28
0.100 Cal Check	0.100	15:58:46
Air Blank	0.000	15:59:16
Subject Test	0.133*	16:02:24
Air Blank	0.000	16:02:56
Air Blank	0.000	16:03:24
Subject Test	XXX**	16:04:04
Air Blank	0.000	16:04:38

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Sof 7
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT *gk*
Serial Number: 80-002265
12/26/2006 16:05:00

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:05:28
Cal Check	0.099	16:05:47
Air Blank	0.000	16:06:16



Operator's Signature

6 of 7
B4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

Handwritten initials

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 20047 AGENCY Tucson Police Department
DATE 12-26-2006 TIME 1544
Intoxilyzer Serial # 80-002265 Location Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 7
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
12/26/2006 15:44:55

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:45:22
Cal Check	0.100	15:45:41
Air Blank	0.000	15:46:10


Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

1 of 1
TS4

Intoxilyzer Serial Number: 80-002265

Initial Service Date: 12-26-2006

CMI, Inc. Certificate of Calibration Date: 11-10-2006

Comment: The regulator leaked when received. A new regulator was sent by CMI, Inc. and installed on 12-26-2006. TS4 28847



AN EMPLOYEE OWNED COMPANY

MPD, INC.

316 East Ninth Street
Owensboro, Kentucky 42303

SHIPPING NOTICE

MARK ALL BILLS OF LADING,
PACKAGES, INVOICES AND
CORRESPONDENCE WITH
THIS NUMBER

No. **PIL121806B**

SHIP TO:

Jennifer Coombs
Tucson Police Department
270 South Stone Avenue
Tucson, AZ 85701-1917

CHECK REASON FOR SHIPMENT:

- 1. Material Returned to Vendor
- 2. Customer Product Returns
- 3. Material @ Vendor Outside Proc.
- 4. Sale of Equipment
- 5. Other _____

VENDOR CODE NO.

ACCOUNT NO.

268

UNIT	QUANTITY	UNIT OF MEASURE	UNIT PRICE	PER	DRAWING NO. AND/OR DESCRIPTION	TOTAL PRICE
1	1		na	1	Gas Regulator Replacement	

CHECK METHOD OF SHIPMENT:

- Parcel Post
- UPS Surface
- UPS Next Day
- UPS 2nd Day
- Truck _____
- Fed. Exp. Std.
- Fed. Exp. Pr. One
- Fed. Exp. Sat. Del.
- Other

CHECK FOR SHIPPING CHARGES:

- Prepaid
- Collect
- Prepay & Add

SHIPPING DEPARTMENT USE ONLY:

Shipped By: _____

Date Shipped: _____

How Shipped: _____

Waybill No. _____

Trans. Chgs.: _____

Pursuant to section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Signature: ✓ Phillip I Lively