



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

Tdy
6.19.07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-7-07 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 5-10-07 and 6-7-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

June 19, 2007

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 19, day of June, 2007.

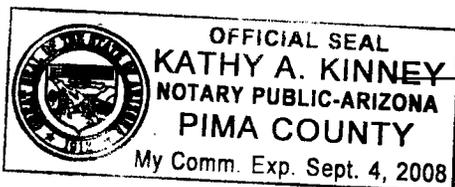

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

~~1 of 8 ca ca~~
1 of 9 ca

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TJ
6-18-07

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 6-7-07 TIME 1446

INTOXILYZER SERIAL # 80-00 2265 LOCATION ODD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS Deficient sample check was not performed as part
of the test sequence. When this was discovered, the
check was performed on 6-12-07.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

~~2 of 8 cu cu~~
2 of 9 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/07/2007 14:41:55

Tsh
6-18-07

Standard Lot#: 635301I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:42:54
Diagnostic Test	Pass	14:43:21
Air Blank	0.000	14:43:48
0.100 Cal Check	0.101	14:44:07
Air Blank	0.000	14:44:37
Subject Test	0.000	14:44:53
Air Blank	ABT*	14:45:11
Air Blank	0.000	14:45:40

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

~~3 of 8 cu cu~~
3 of 9 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/07/2007 14:47:36

TS4
6-18-07

Standard Lot#: 635301I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:48:33
Diagnostic Test	Pass	14:49:01
Air Blank	IPS*	14:49:23
Air Blank	0.000	14:49:51

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

~~4 of 8 cu cu~~
4 of 9 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/07/2007 14:49:56

TPH
6.18.07

Standard Lot#: 635301I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:50:56
Diagnostic Test	Pass	14:51:23
Air Blank	0.000	14:51:51
0.100 Cal Check	RFI*	14:52:03
Air Blank	0.000	14:52:33

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

~~5068 cu cu~~
5069 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/07/2007 14:52:54

Standard Lot#: 635301I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

TS4
6-18-07

Test	g/210L	Time
Air Blank	0.000	14:53:49
Diagnostic Test	Pass	14:54:17
Air Blank	0.000	14:54:44
0.100 Cal Check	0.100	14:55:03
Air Blank	0.000	14:55:33
Subject Test	XXX*	14:57:24
Air Blank	0.000	14:57:59

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
06/07/2007 14:58:32

~~6 of 8 in~~
6 of 9 in

Standard Lot#: 635301I
Last Changed By: T. GALLEGOS #28847

TS4
6-18-07

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:59:00
Cal Check	0.100	14:59:19
Air Blank	0.000	14:59:48

Clare Conley
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

~~7 of 8~~ cu
7 of 9 cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4 6.807

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 6-7-07 TIME 1446
Intoxilyzer Serial # 80-002265 Location 0DD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
06/07/2007 14:46:09

Standard Lot#: 635301I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:46:37
Cal Check	0.100	14:46:56
Air Blank	0.000	14:47:25

Clare Conley

Operator's Signature

8 of 800
or

8 of 9 cu

TJ
6-18-07

9 of 9 cu

TSh
6-18-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/12/2007 13:36:38

Standard Lot#: 635301I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:37:34
Diagnostic Test	Pass	13:38:02
Air Blank	0.000	13:38:29
0.100 Cal Check	0.100	13:38:48
Air Blank	0.000	13:39:17
Subject Test	0.083*	13:42:26
Air Blank	0.000	13:42:57
Air Blank	0.000	13:43:25
Subject Test	ABT**	13:43:35
Air Blank	0.000	13:44:04

*Deficient Sample
**Sequence Aborted

Not a Successfully Completed Test Sequence