



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

CU
11-23-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10.25.2007 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 9.26.2007 and 10.25.2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .


Terry S. Gallegos

Terry S. Gallegos

11.21.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 21, day of November, 2007.

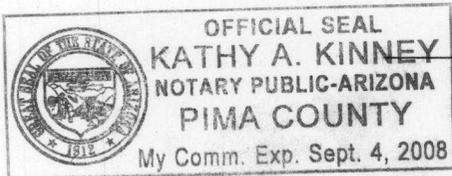

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

W
11.23.07

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

10 of 7
B4

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 20847 AGENCY Tucson Police Department

DATE 10-25-2007 TIME 1054

INTOXILYZER SERIAL # 80-002265 LOCATION ODD

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓✓
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

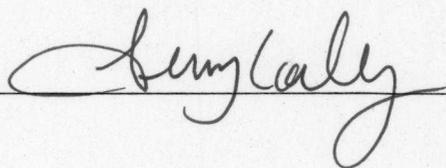SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
10/25/2007 10:56:06

cu

2 of 7
BY

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:56:55
Diagnostic Test	Pass	10:57:23
Air Blank	0.000	10:57:50
0.100 Cal Check	0.100	10:58:09
Air Blank	0.000	10:58:38
Subject Test	0.000	10:59:21
Air Blank	0.000	10:59:48
Five Minute Wait	ABT*	
Air Blank	0.000	11:00:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
10/25/2007

11:00:40

CU
30 of 7
754

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:01:30
Diagnostic Test	Pass	11:01:57
Air Blank	0.000	11:02:24
0.100 Cal Check	RFI*	11:02:37
Air Blank	IPS**	11:02:46

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
10/25/2007

11:02:57

cu
4 of 7
B4

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:03:48
Diagnostic Test	Pass	11:04:15
Air Blank	0.000	11:04:43
0.100 Cal Check	0.100	11:05:02
Air Blank	0.000	11:05:31
Subject Test	0.253*	11:08:39
Air Blank	0.000	11:09:11
Air Blank	0.000	11:09:39
Subject Test	XXX**	11:10:03
Air Blank	0.000	11:10:37

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

10/25/2007

11:10:55


Standard Lot#: 719801I

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	11:11:22
Cal Check	0.100	11:11:41
Air Blank	0.000	11:12:11

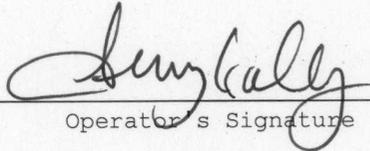

Operator's Signature

u
5 of 7
B4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)


ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

60f 711.23-07
TS4

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 20047 AGENCY Tucson Police Department
DATE 10-25-2007 TIME 1054
Intoxilyzer Serial # 80-002265 Location ODD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

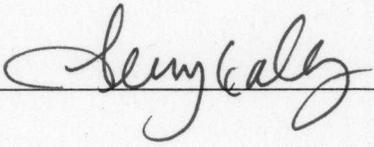SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

10/25/2007

10:54:14

cu

7077
TSy


Standard Lot#: 719801I

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	10:54:42
Cal Check	0.100	10:55:00
Air Blank	0.000	10:55:30

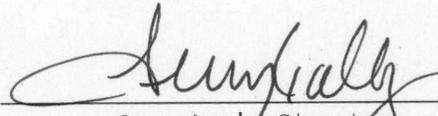

Operator's Signature