



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY
12-21-07

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/14/07 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 11/21/07 and 12/14/07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah L. Cione

Sarah L. Cione

12/19/07

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19, day of December, 2007.

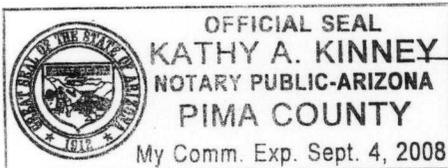

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 SC
134

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 12/14/07 TIME 15:44
INTOXILYZER SERIAL # 80-002265 LOCATION Team 5

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
12/14/2007

15:45:58

2 of 7sc

TS4

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:46:51
Diagnostic Test	Pass	15:47:18
Air Blank	0.000	15:47:46
0.100 Cal Check	0.100	15:48:04
Air Blank	0.000	15:48:34
Subject Test	0.000	15:48:59
Air Blank	0.000	15:49:26
Five Minute Wait	ABT*	
Air Blank	0.000	15:50:01

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
12/14/2007 15:50:07

B4

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	15:50:43
Air Blank	IPS**	15:50:53

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
12/14/2007 15:51:04

4 of 7 SC

BY

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:51:53
Diagnostic Test	Pass	15:52:20
Air Blank	0.000	15:52:47
0.100 Cal Check	0.101	15:53:06
Air Blank	0.000	15:53:36
Subject Test	0.170*	15:56:44
Air Blank	0.000	15:57:16
Air Blank	0.000	15:57:44
Subject Test	XXX**	15:58:11
Air Blank	0.000	15:58:44

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
12/14/2007 15:59:04

5 of 7 sc
134

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:59:31
Cal Check	0.100	15:59:50
Air Blank	0.000	16:00:19

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 SC

154

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 12/14/07 TIME 15:44
Intoxilyzer Serial # 80-002265 Location Team 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
12/14/2007 15:44:27

7 of 7sc
154

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:44:55
Cal Check	0.100	15:45:14
Air Blank	0.000	15:45:43

Sarah Cione

Operator's Signature