



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

Tby  
1-29-2008

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/28/08 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 1/11/08 and 1/28/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

1/29/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 29, day of January, 2008.

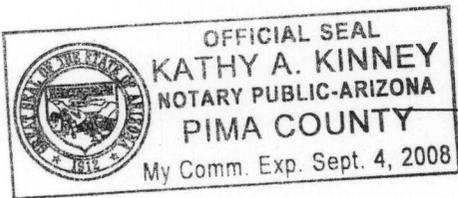

Kathy A Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

By  
1-29-2008

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/28/08 TIME 14:36

INTOXILYZER SERIAL # 80-002265 LOCATION Team 5

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS (26227)  
Officer George Eppley reported a Diagnostic Failure - Internal Printer  
Test. Performed SQAP to check Diagnostic Failure

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 75C

104

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
01/28/2008 14:37:33

Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:38:22
Diagnostic Test	Pass	14:38:49
Air Blank	0.000	14:39:16
0.100 Cal Check	0.100	14:39:35
Air Blank	0.000	14:40:05
Subject Test	0.000	14:40:29
Air Blank	0.000	14:40:56
Five Minute Wait	ABT*	
Air Blank	0.000	14:41:32

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 75C

BY

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
01/28/2008 14:42:29

Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:43:22
Diagnostic Test	Pass	14:43:50
Air Blank	0.000	14:44:17
0.100 Cal Check	0.100	14:44:36
Air Blank	0.000	14:45:05
Subject Test	0.261*	14:48:13
Air Blank	0.000	14:48:46
Air Blank	0.000	14:49:13
Subject Test	XXX**	14:49:35
Air Blank	0.000	14:50:08

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

4 of 7 SC  
B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
01/28/2008 14:50:29

---

Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:51:27
Diagnostic Test	Pass	14:51:55
Air Blank	IPS*	14:52:20
Air Blank	RFI**	14:52:29

\*Improper Sample  
\*\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
01/28/2008 14:52:48

5 of 7 sc

PH

-----  
Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:53:16
Cal Check	0.099	14:53:35
Air Blank	0.000	14:54:04

*Sarah Cione*

Operator's Signature

6 of 7 SC  
B4

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 1/28/08 TIME 14:36  
Intoxilyzer Serial # 80-002265 Location Team 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
01/28/2008 14:36:05

7 of 7sc  
B4

-----  
Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:36:32
Cal Check	0.100	14:36:51
Air Blank	0.000	14:37:20

*Sarah Cione*

Operator's Signature