



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
11.3.08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/17/08 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 9/19/08 and 10/17/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/20/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20th, day of October, 2008.

A.M. Hernandez

Notary Public

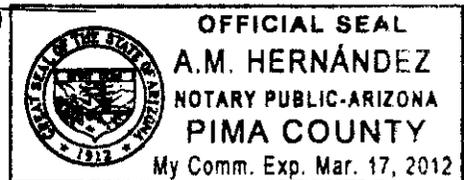

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 75C

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

By
11.3.08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/17/08 TIME 14:19

INTOXILYZER SERIAL # 80-002265 LOCATION Team 5

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7 SC

BS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
10/17/2008 14:20:53

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:21:46
Diagnostic Test	Pass	14:22:13
Air Blank	0.000	14:22:40
0.100 Cal Check	0.101	14:22:59
Air Blank	0.000	14:23:29
Subject Test	0.000	14:23:45
Air Blank	0.000	14:24:12
Five Minute Wait	ABT*	
Air Blank	0.000	14:24:46

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 75c

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
10/17/2008 14:24:51

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:25:19
Air Blank	IPS**	14:25:27

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7
SC

BY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
10/17/2008 14:25:37

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:26:22
Diagnostic Test	Pass	14:26:49
Air Blank	0.000	14:27:16
0.100 Cal Check	0.101	14:27:35
Air Blank	0.000	14:28:05
Subject Test	0.060*	14:31:13
Air Blank	0.000	14:31:45
Air Blank	0.000	14:32:12
Subject Test	XXX**	14:32:39
Air Blank	0.000	14:33:12

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
10/17/2008 14:33:27

5 of 7_{sc}
TS4

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:33:54
Cal Check	0.100	14:34:13
Air Blank	0.000	14:34:42



Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
11-3-08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/17/08 TIME 14:19
Intoxilyzer Serial # 80-002265 Location Team 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
10/17/2008 14:19:07

7 of 75c

TS4

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:19:34
Cal Check	0.100	14:19:53
Air Blank	0.000	14:20:22

Sarah Cione
Operator's Signature