



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
2-4-2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/9/09 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 12/12/08 and 1/9/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/12/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 12th, day of January, 2009.

A.M. Hernandez

Notary Public

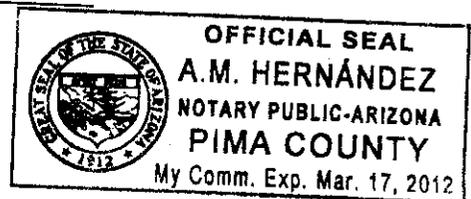

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 9 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
2.4.2009

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 1/9/09 TIME 14:09
INTOXILYZER SERIAL # 80-002265 LOCATION Team 5-0DD

- (✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 9 sc

BY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
01/09/2009

14:12:24

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:13:16
Diagnostic Test	Pass	14:13:43
Air Blank	0.000	14:14:11
0.100 Cal Check	0.101	14:14:29
Air Blank	0.000	14:14:59
Subject Test	0.000	14:15:15
Air Blank	0.000	14:15:42
Five Minute Wait	ABT*	
Air Blank	0.000	14:16:36

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 9 sc

BH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
01/09/2009

14:16:40

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:17:33
Air Blank	IPS**	14:17:43

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 9 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
01/09/2009

14:17:47

BY

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:18:32
Diagnostic Test	Pass	14:18:59
Air Blank	0.000	14:19:26
0.100 Cal Check	0.101	14:19:45
Air Blank	0.000	14:20:15
Subject Test	0.027*	14:23:23
Air Blank	0.000	14:23:54
Air Blank	0.000	14:24:22
Subject Test	XXX**	14:24:40
Air Blank	0.000	14:25:12

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 9 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
01/09/2009

14:26:03

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

134

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:26:31
Cal Check	0.100	14:26:50
Air Blank	0.000	14:27:19

Sarah Cione
Operator's Signature

EXHIBIT G-2

6 of 9 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
2.4.2009

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/9/09 TIME 14:09
Intoxilyzer Serial # 80-002265 Location Team 5-0DD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 9 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
01/09/2009 14:09:50

TS4

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:10:18
Cal Check	0.101	14:10:37
Air Blank	0.000	14:11:06

Sarah Cione
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 1/9/09
QAS: Sarah Cione
Intoxilyzer Serial Number: 80-002265
New Lot Number: 801002 I
New Expiration Date: 1/15/2010

TSY
2-4-09

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

9 of 9 SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
01/09/2009

14:31:11

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:31:39
Cal Check	0.101	14:31:57
Air Blank	0.000	14:32:27

Sarah Cione

Operator's Signature